


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ritodrine

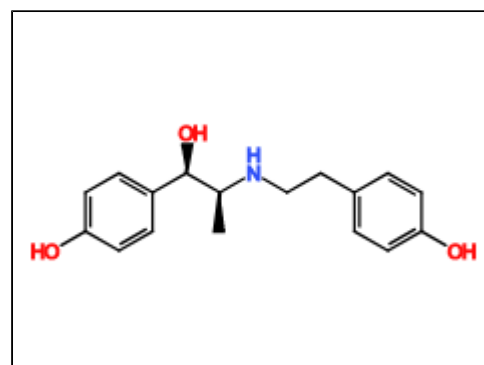

CHROMATOGRAM

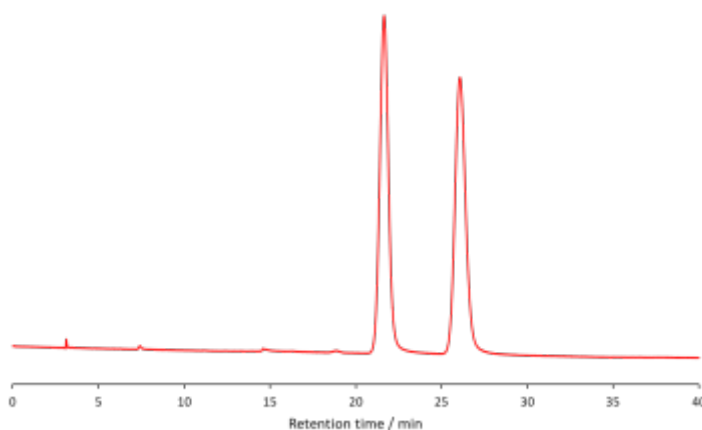

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (μL)		Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	21.6	5.93	1 - 2	4.30
2	26.1	7.35		