


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Afatinib

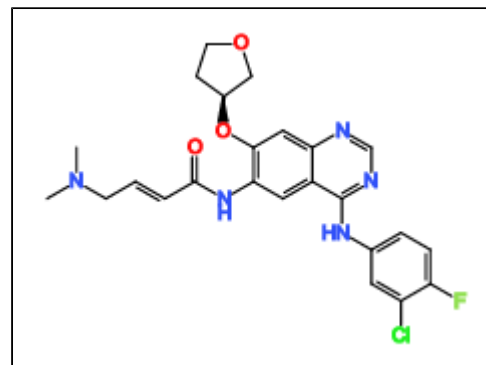

CHROMATOGRAM

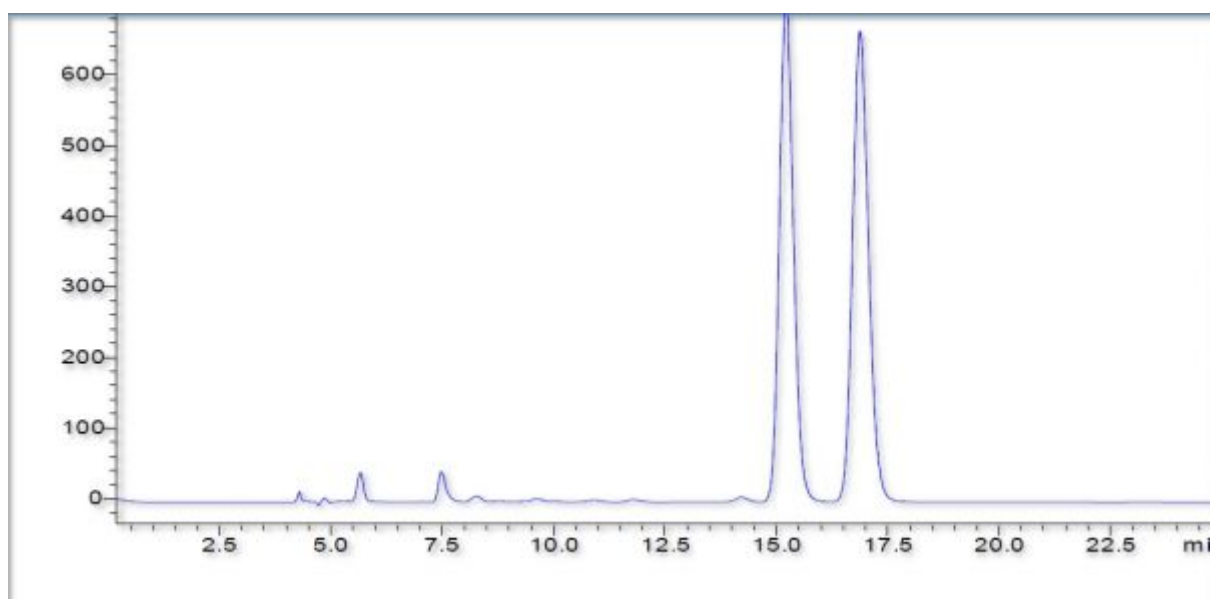

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IE 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether / methanol / diethylamine = 20 / 80 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	15
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 340 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	15.2		1 - 2	2.50
2	17.0			