


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Alvimopan

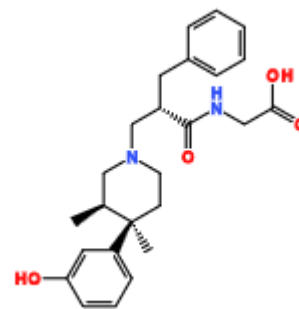

CHROMATOGRAM

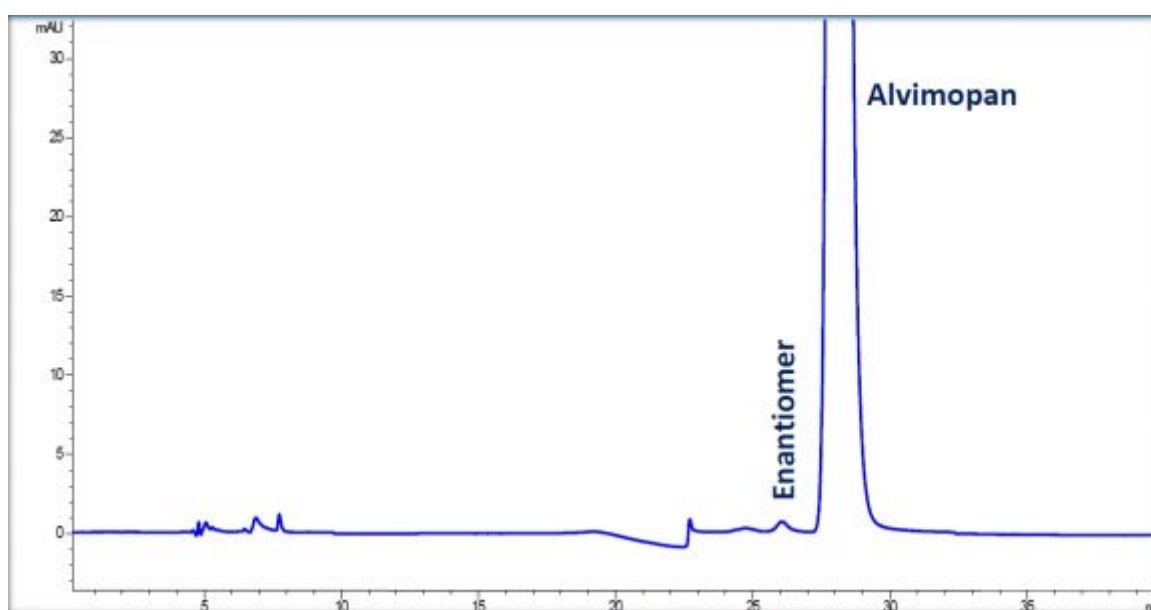

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	4 x 250 mm (3 µm)		
Mobile Phase Composition	acetonitrile / tetrahydrofuran / water / formic acid / diethylamine = 490 / 490 / 20 / 1.9 / 2.6		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 273 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	26.1		1 - 2	2.63
2	28.2			