


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Batefenterol

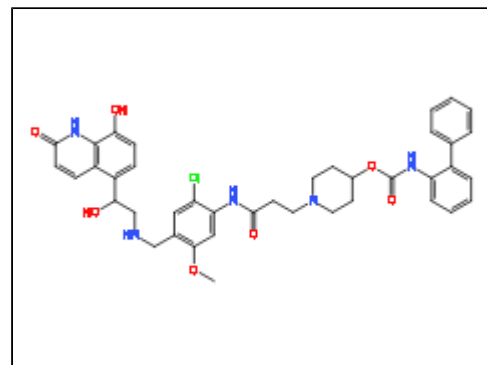

CHROMATOGRAM

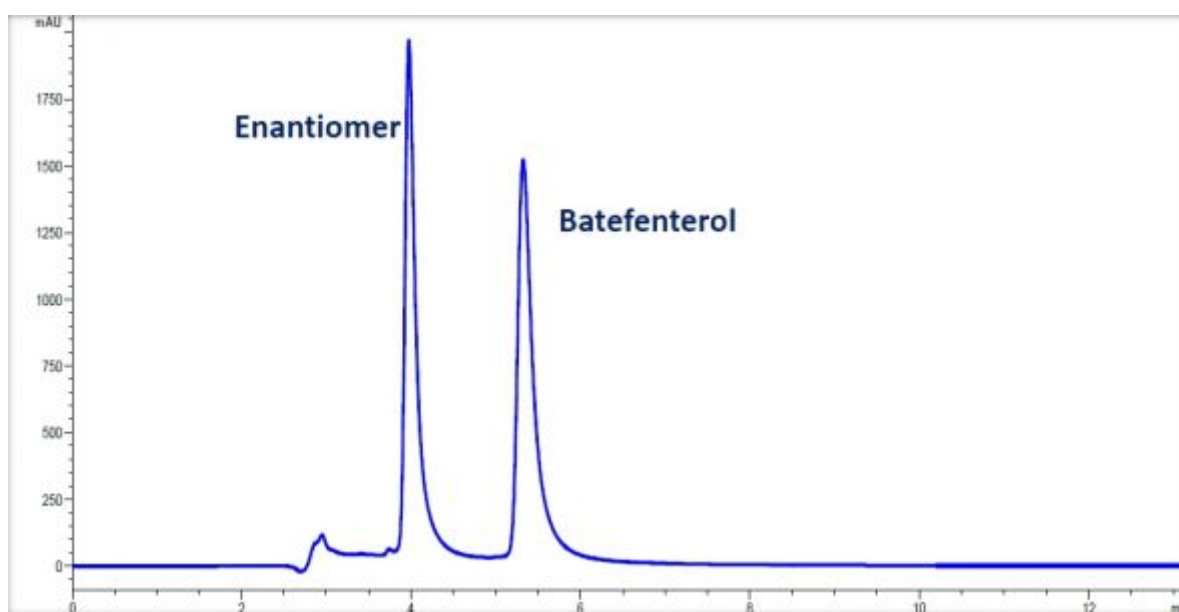

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(-) 3 x 250 mm (3 µm)		
Mobile Phase Composition	methanol / tetrahydrofuran / water / formic acid / diethylamine = 500 / 500 / 10 / 1.9 / 5.2		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.0		1 - 2	5.00
2	5.3			