


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Carfilzomib

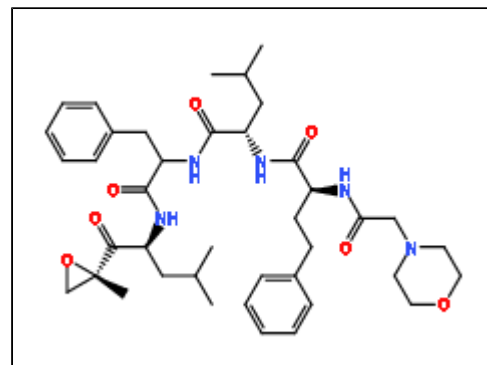

CHROMATOGRAM

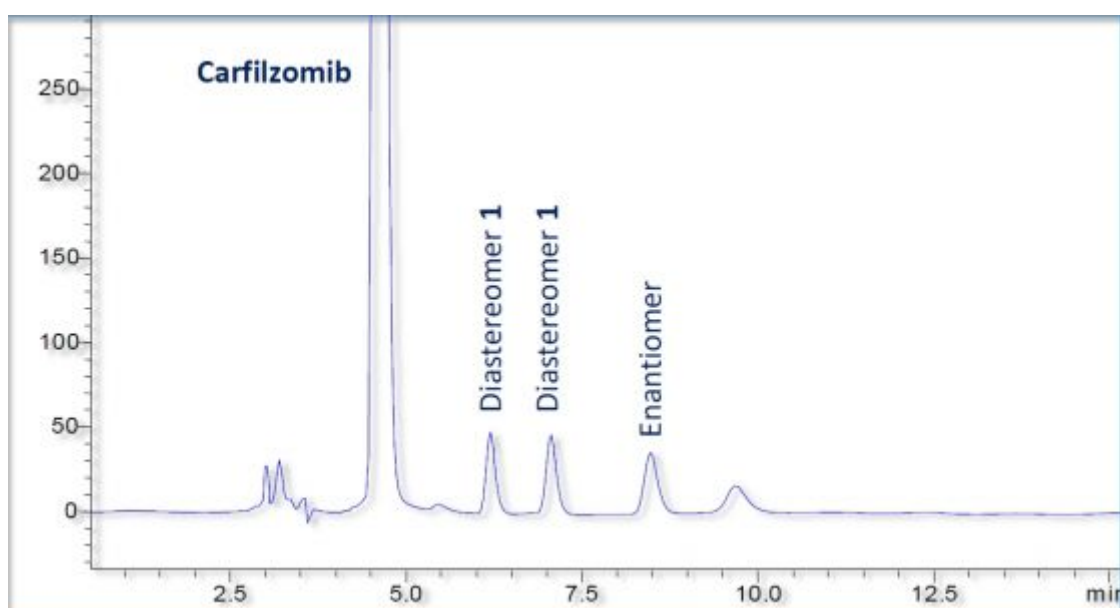

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / methyl tert-butyl ether / methanol = 50 / 40 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	30
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.6		1 - 2	5.30
2	6.2		2 - 3	2.80
3	7.0		3 - 4	3.90
4	8.5			