


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Droxidopa

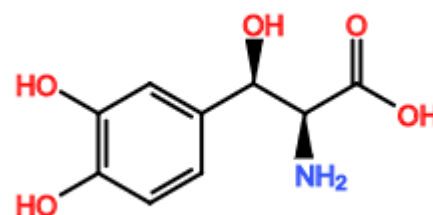

CHROMATOGRAM

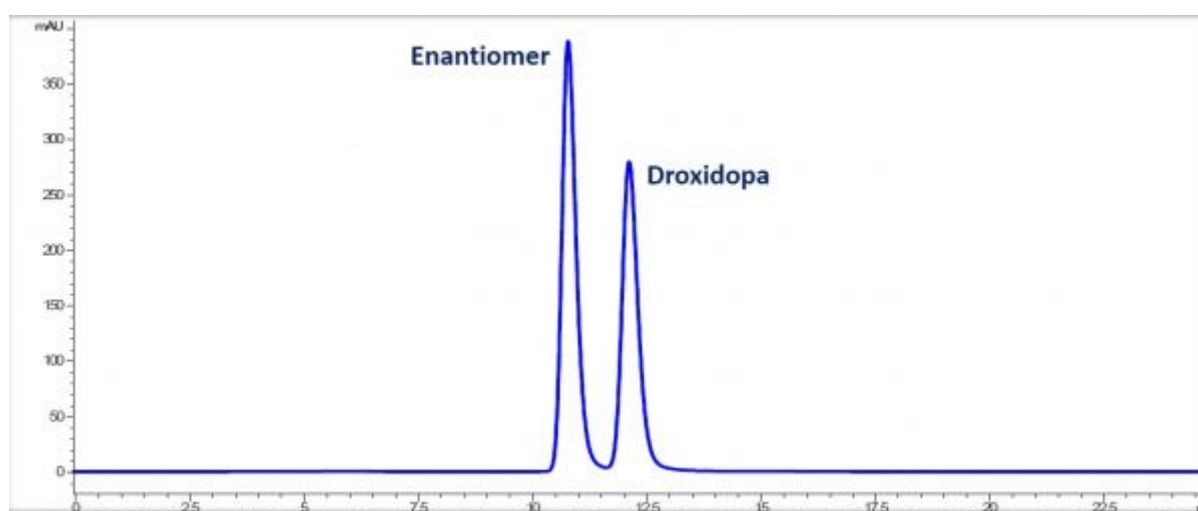

ANALYTICAL REPORT

Column Trade Name	CROWNPAK [®] CR(+)		
Column Size (Particle Size)	4 x 150 mm (5 μm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 = 100		
Chromatographic Mode	Reverse Phase		

Flow rate (mL/min)	0.2	Temperature (°C)	10
Injection Volume (μL)		Sample Concentration	
Detection	UV-VIS 280 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.8		1 - 2	2.20
2	12.1			