


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Brinzolamide

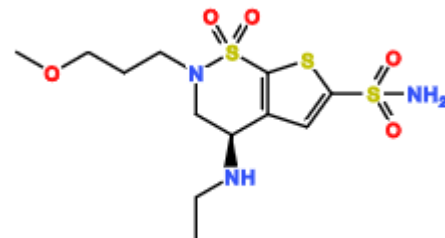

## CHROMATOGRAM

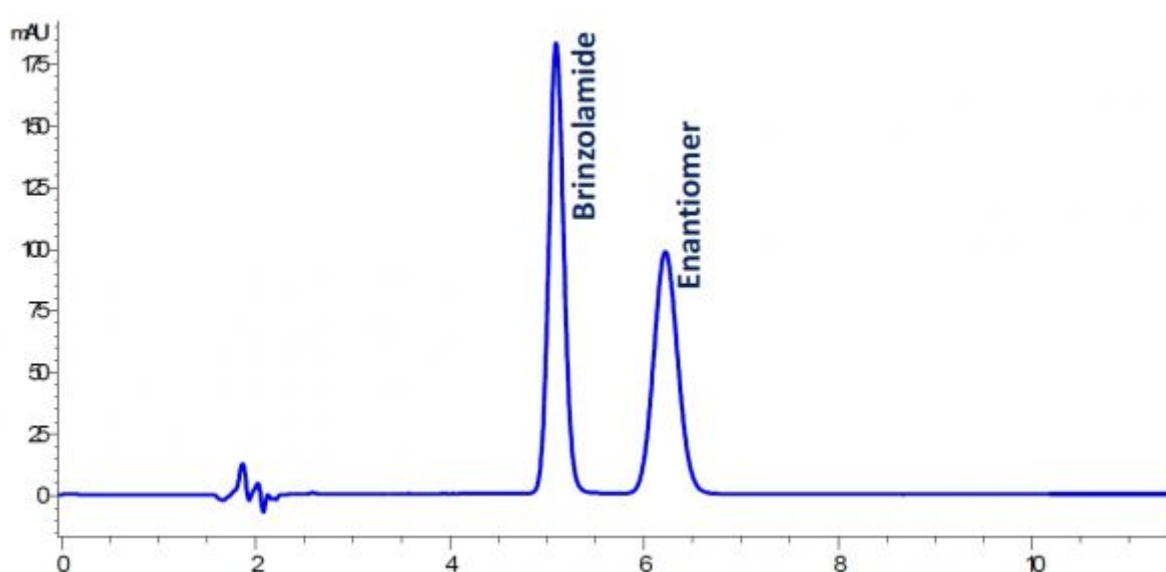

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IG-3 4.6 x 150 mm (3 µm)		
Mobile Phase Composition	10mM ammonium bicarbonate aq. / acetonitrile = 70 / 30		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.1		1 - 2	2.99
2	6.2			