


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Brompheniramine maleate

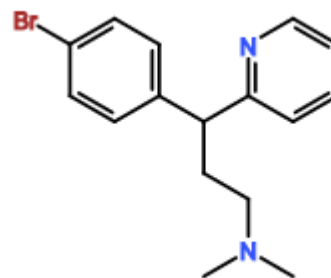

CHROMATOGRAM

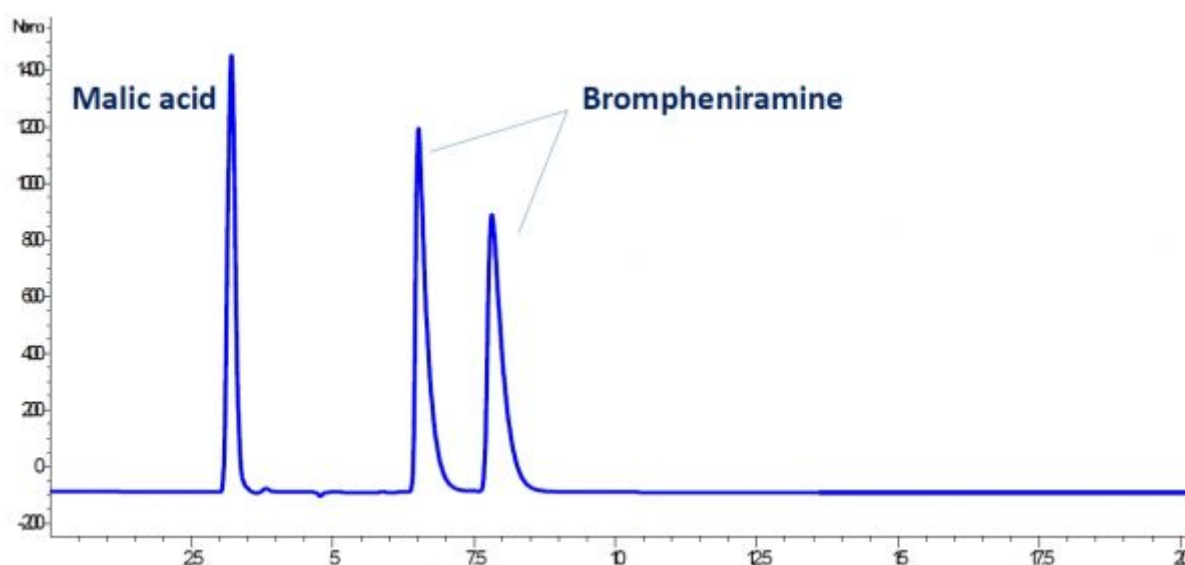

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IG-3 4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / ethanol / ammonium bicarbonate (weight) = 50 / 950 / 0.4		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	15
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 222 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.5		1 - 2	3.40
2	7.8			