


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Safinamide

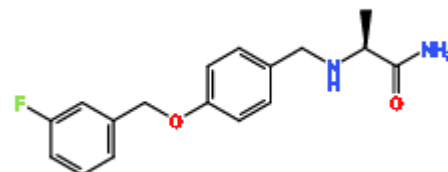

CHROMATOGRAM

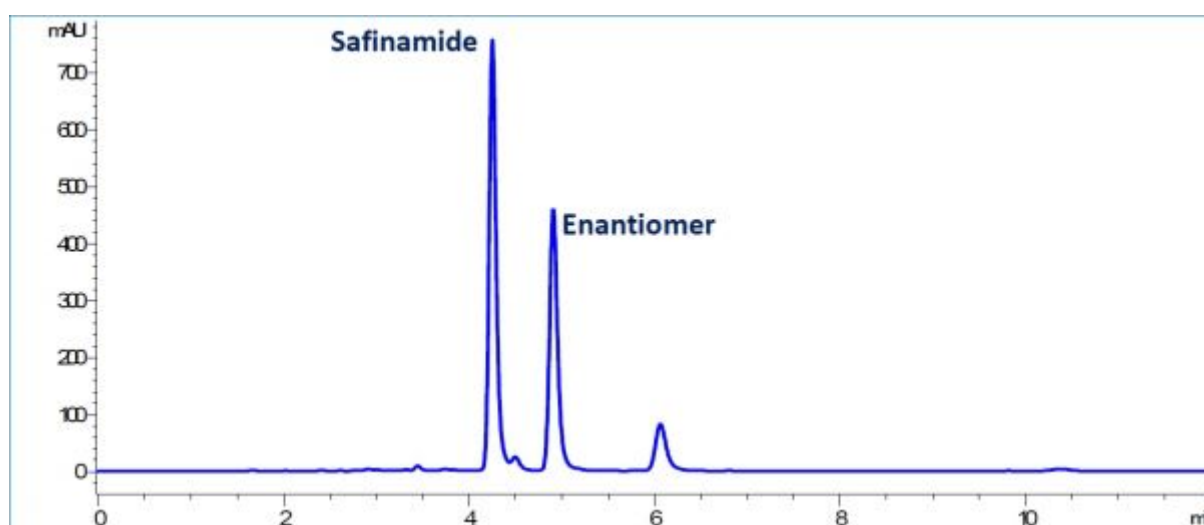

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile / ammonium bicarbonate (weight) = 50 / 950 / 0.4		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	15
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 226 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.3		1 - 2	4.50
2	4.9			