


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Colchicine

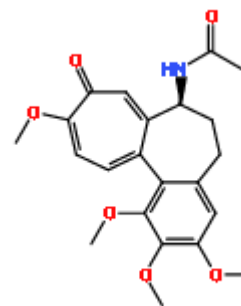

CHROMATOGRAM

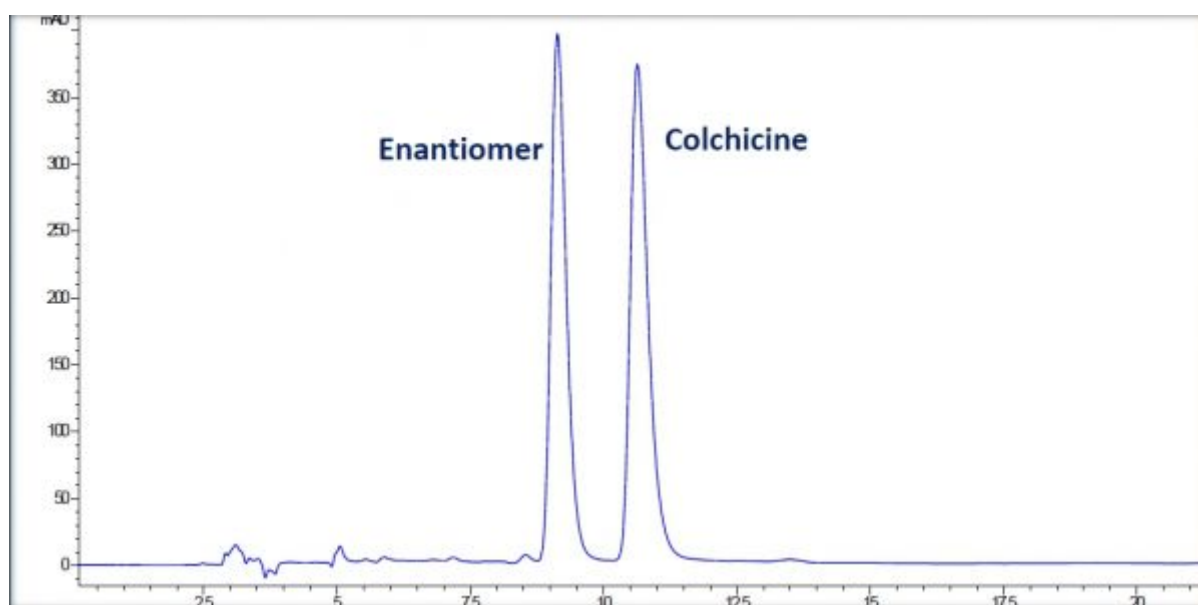

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IG-3 4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 70 / 20 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 245 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.1		1 - 2	2.64
2	10.6			