


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Crizotinib

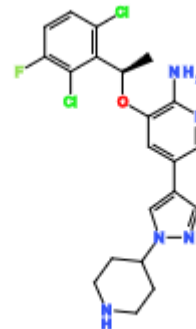

CHROMATOGRAM

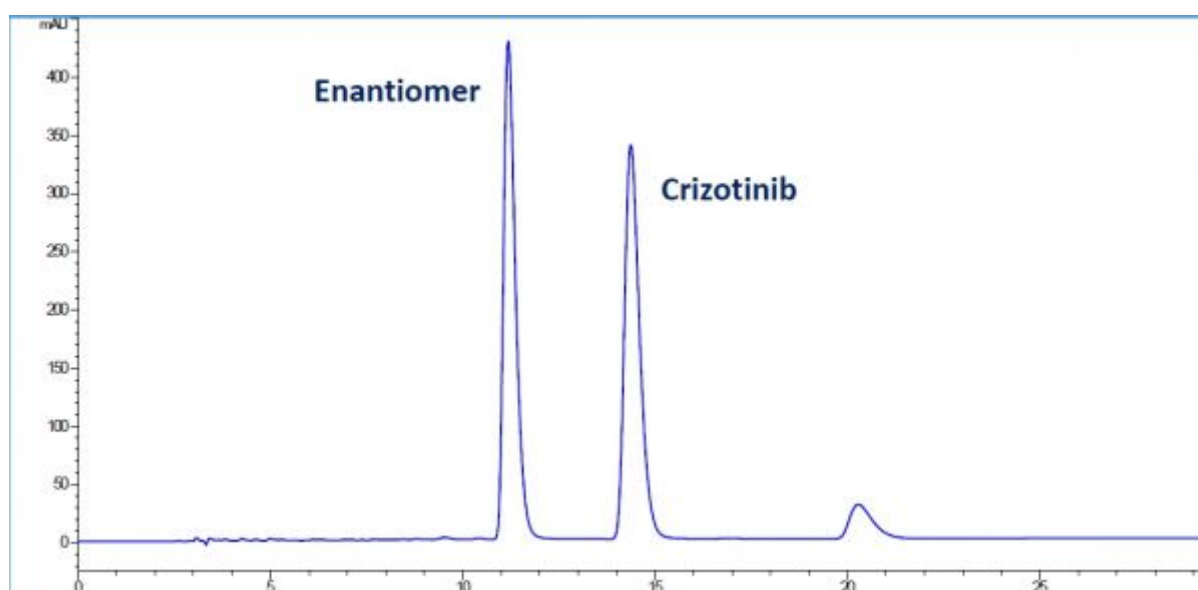

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / methanol / isopropylamine = 50 / 30 / 20 / 0.3		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 265 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.2		1 - 2	4.90
2	14.4			