


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Darunavir

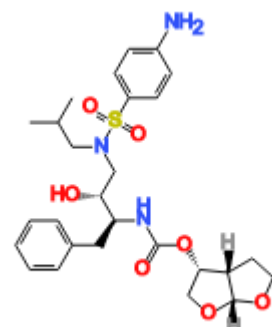

## CHROMATOGRAM

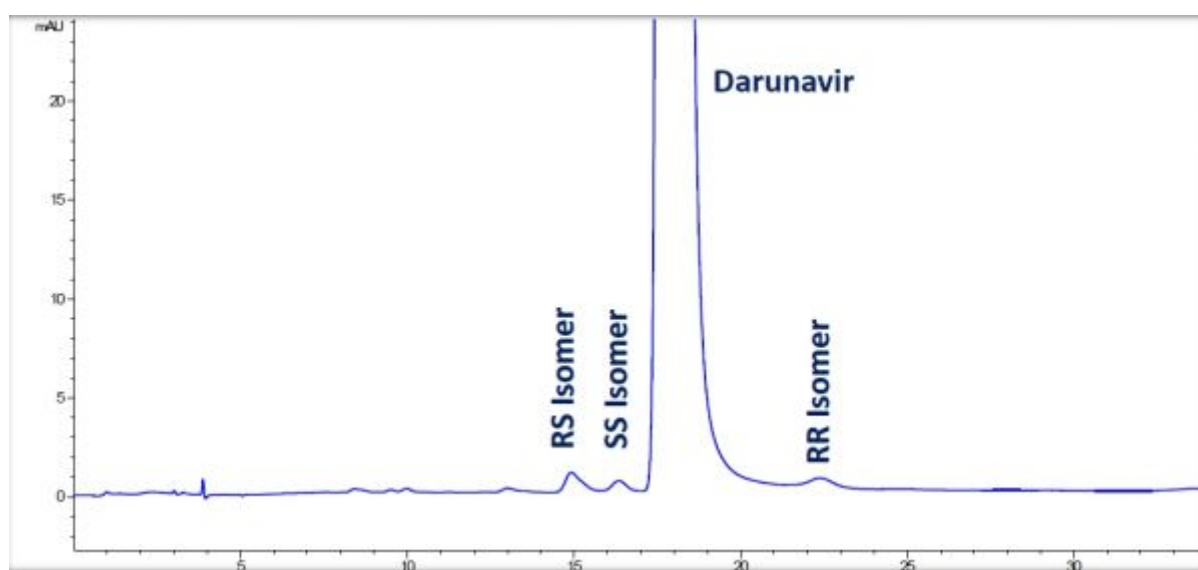

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 265 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	14.9		1 - 2	1.70
2	16.4		2 - 3	1.80
3	17.8		3 - 4	4.30
4	22.4			