


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

5,5-Diphenyl-4-methyl-2-oxazolidinone

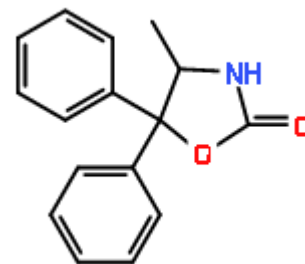

CHROMATOGRAM

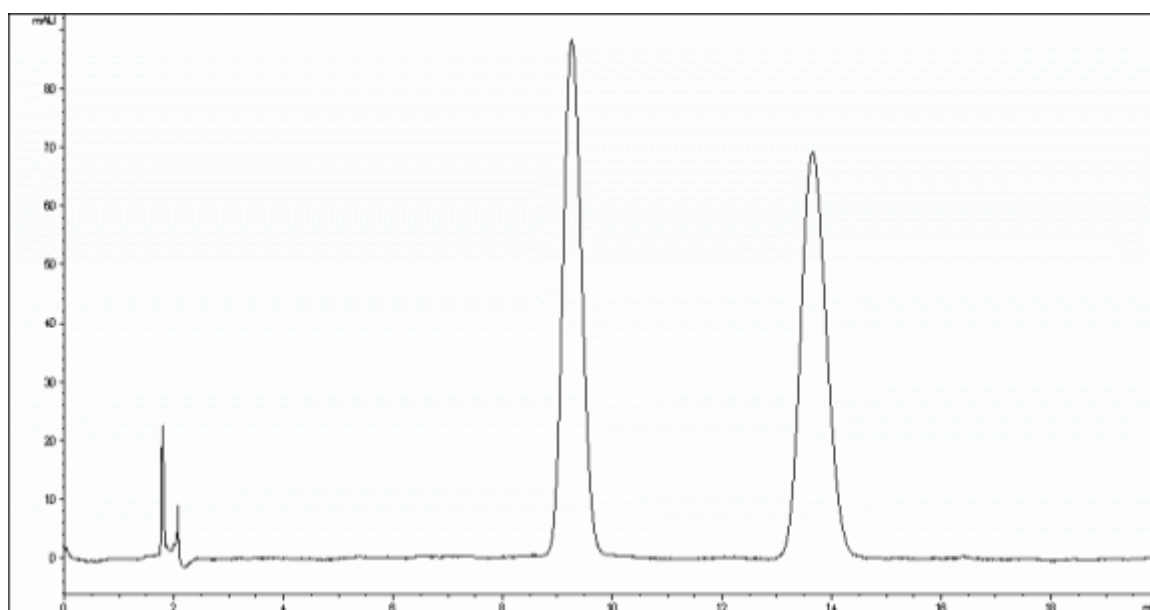

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	methyl tert-butyl ether / ethanol = 98 / 2		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.3	4.42	1 - 2	5.86
2	13.7	6.99		