


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Eletriptan

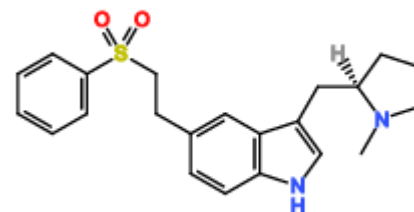

## CHROMATOGRAM

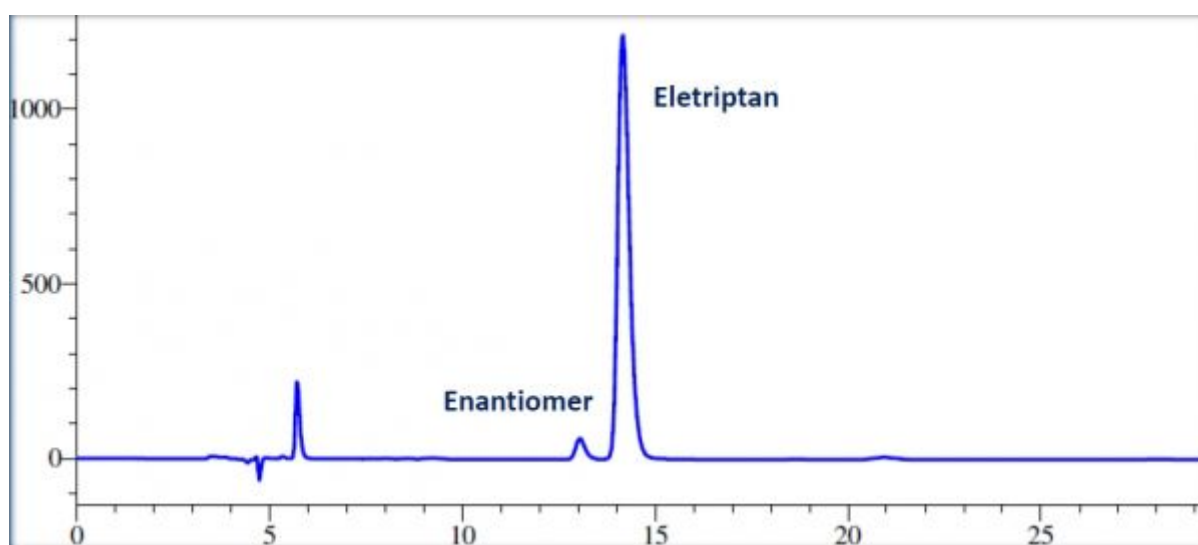

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 70 / 10 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.8	Temperature (°C)	20
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	13.0		1 - 2	2.00
2	14.1			