


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Hesperidin

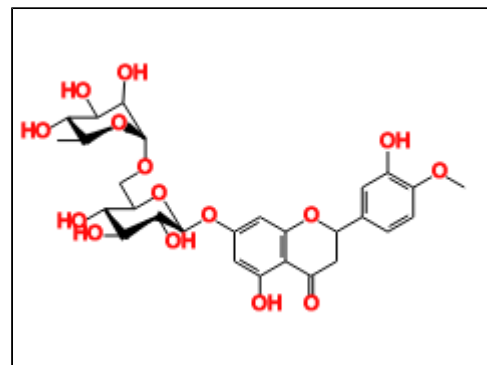

CHROMATOGRAM

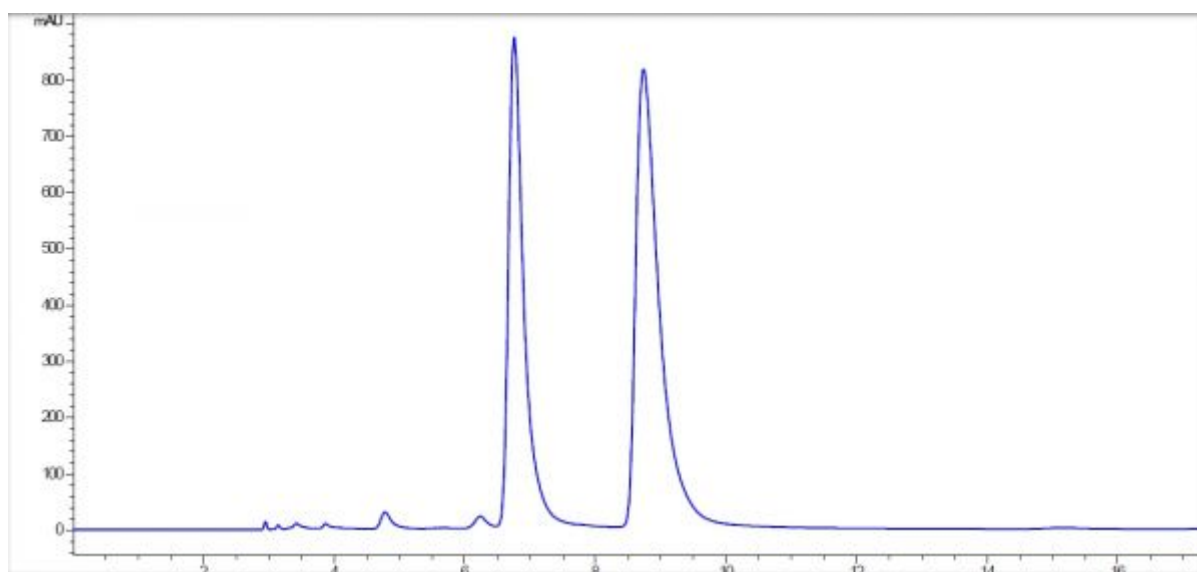

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether / methanol = 85 / 15		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 280 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.8		1 - 2	3.90
2	8.7			