


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Lenalidomide

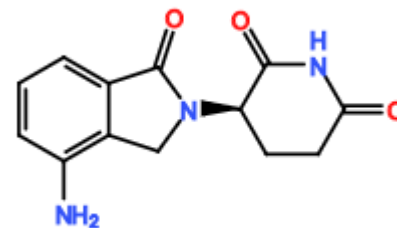

CHROMATOGRAM

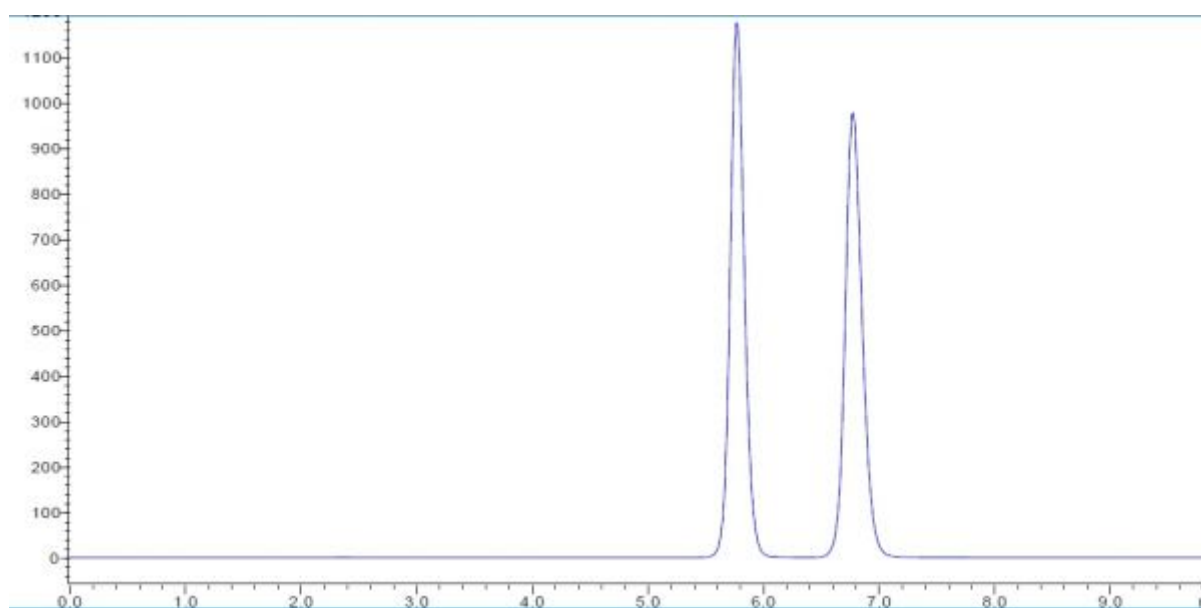

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile / ammonium bicarbonate (weight) = 50 / 950 / 0.79		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	40
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 245 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.8		1 - 2	3.90
2	6.8			