


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Linagliptin

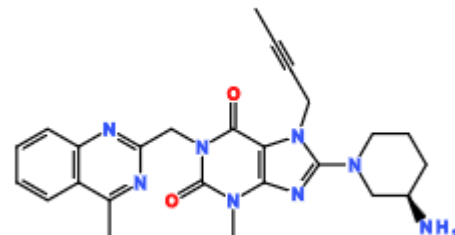

CHROMATOGRAM

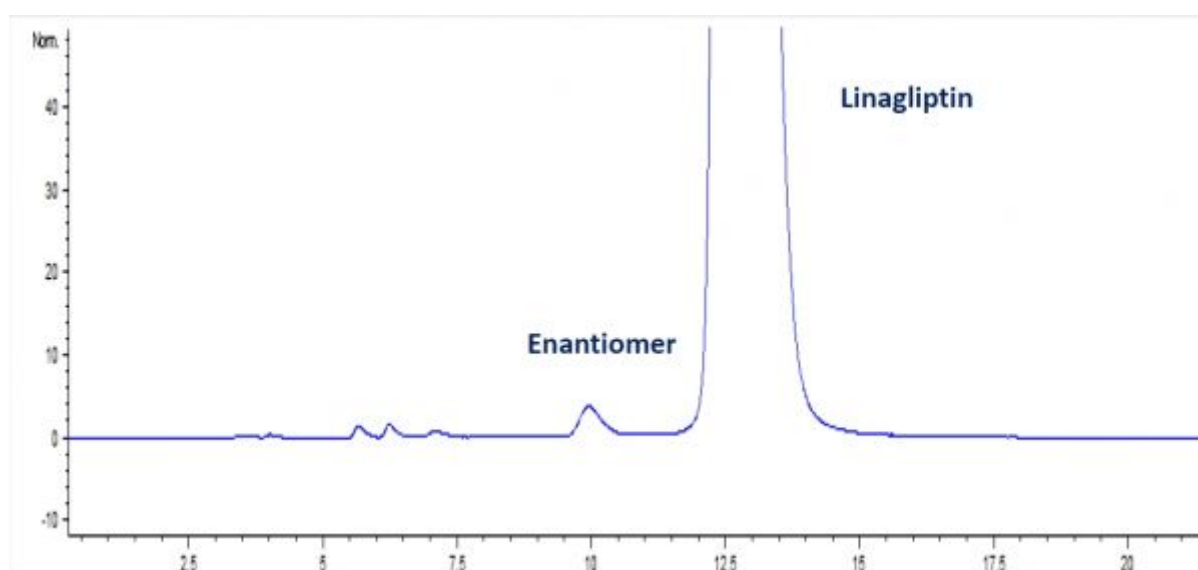

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	ethanol / ethanolamine = 100 / 0.3		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 295 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.9		1 - 2	2.90
2	12.6			