


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Lorcaserin

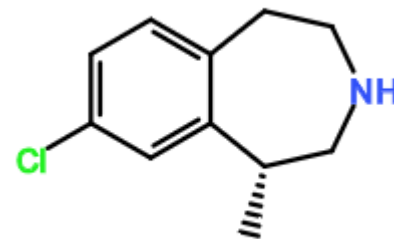

CHROMATOGRAM

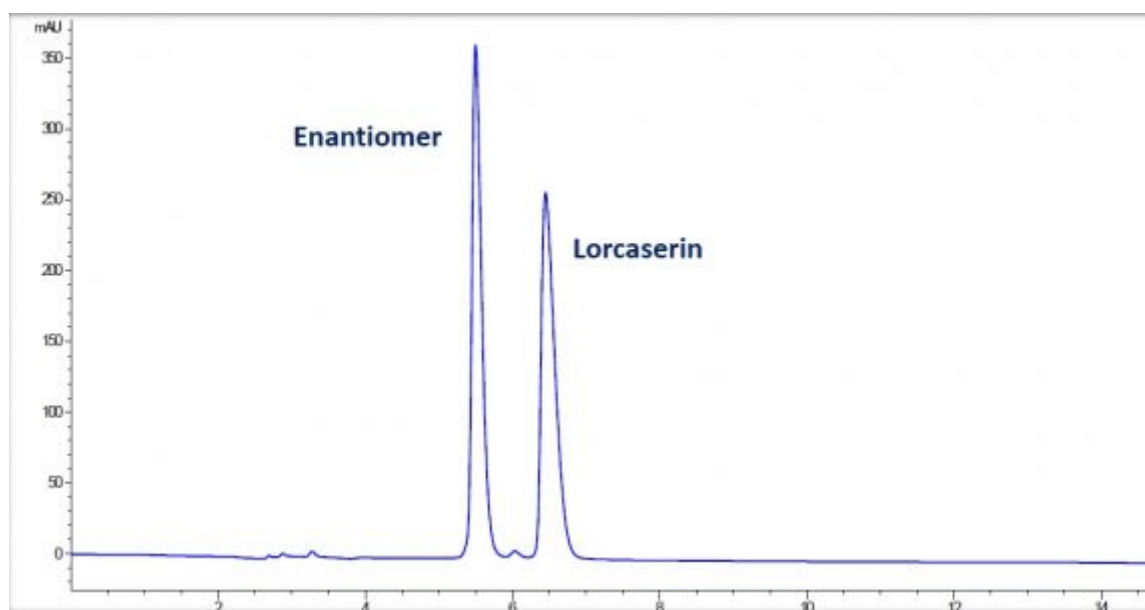

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methanol / ethylenediamine = 100 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.5		1 - 2	3.30
2	6.5			