


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Metaraminol

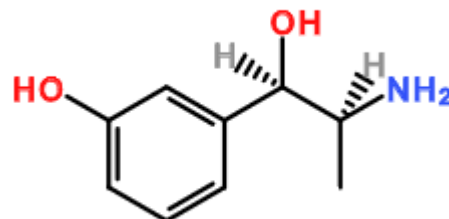

CHROMATOGRAM

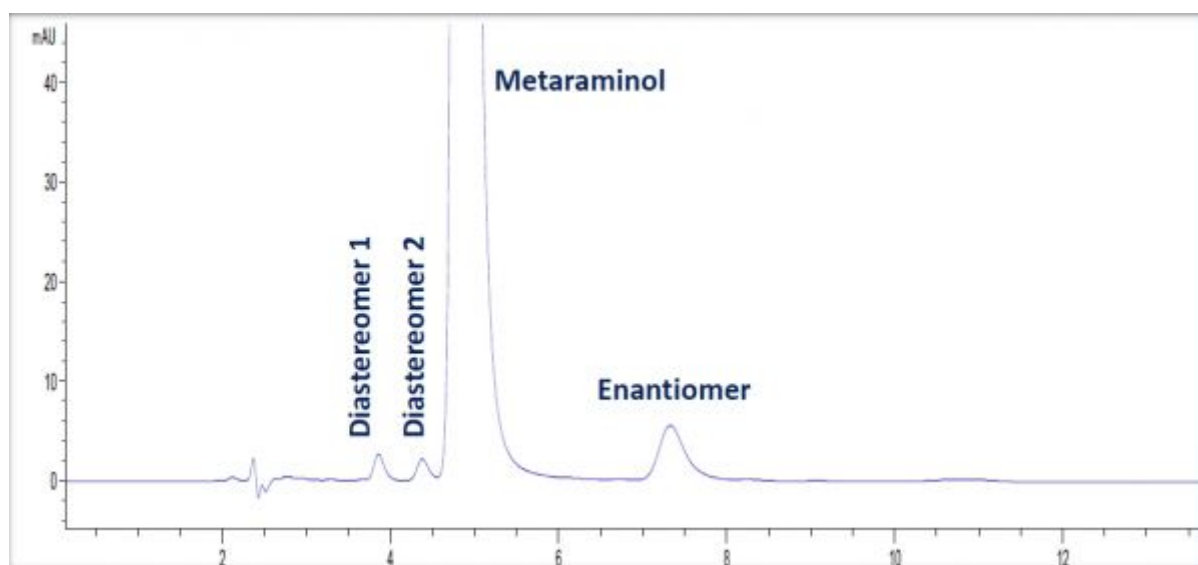

ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.4	Temperature (°C)	10
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.9		1 - 2	2.20
2	4.4		2 - 3	1.60
3	4.8		3 - 4	5.60
4	7.3			