


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Pemetrexed disodium

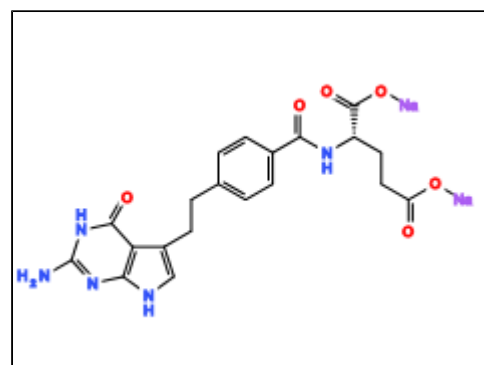

## CHROMATOGRAM

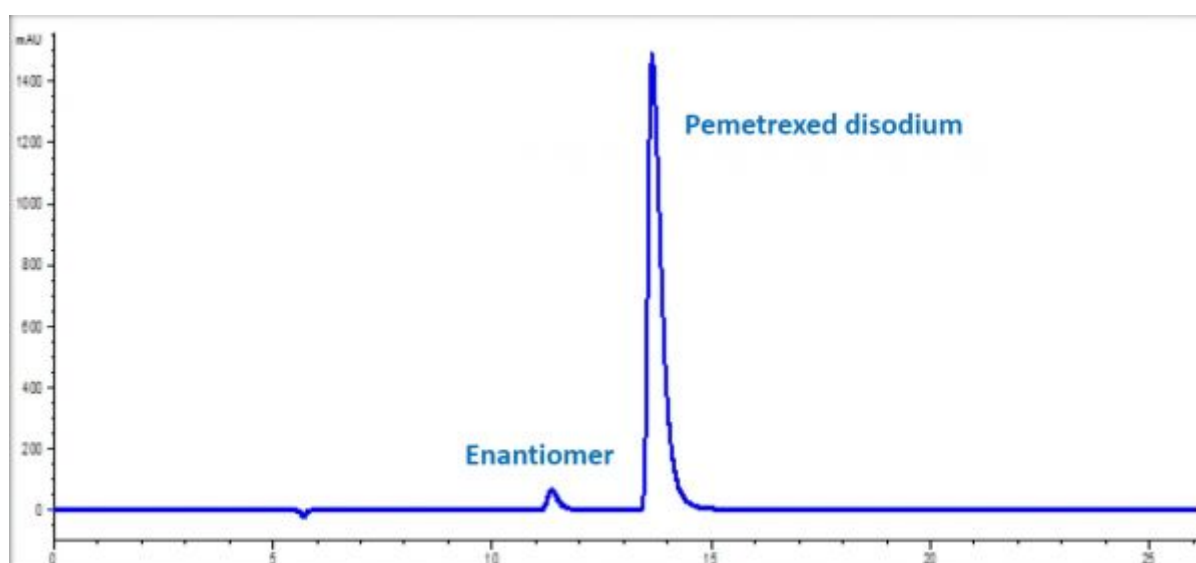

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 4 x 250 mm (3 µm)		
Mobile Phase Composition	methanol / water / formic acid / diethylamine = 980 / 20 / 1.9 / 2.6		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.3		1 - 2	4.73
2	13.6			