


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Posaconazole

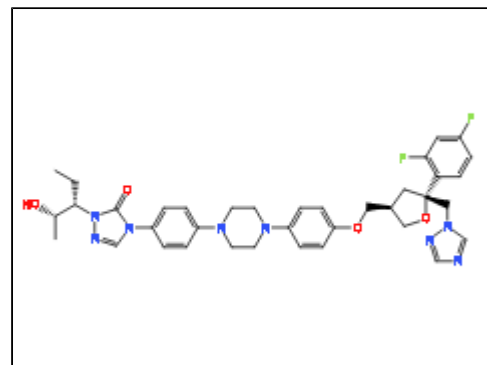

## CHROMATOGRAM

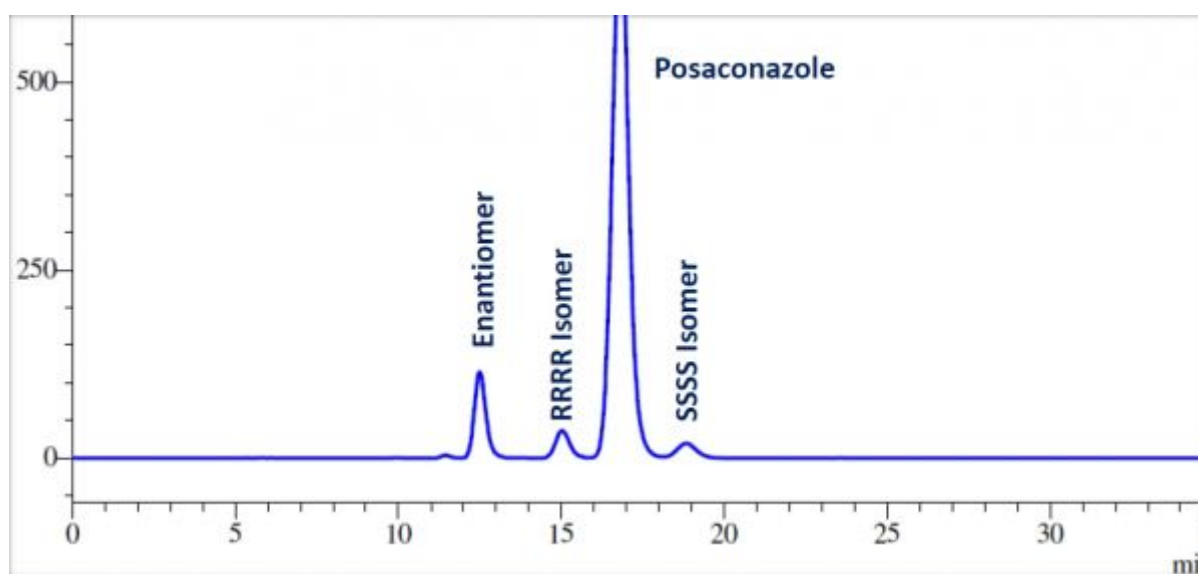

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	dichloromethane / 2-propanol / diethylamine = 50 / 50 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.5		1 - 2	3.40
2	15.0		2 - 3	2.00
3	16.8		3 - 4	1.90
4	18.8			