


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Pregabalin

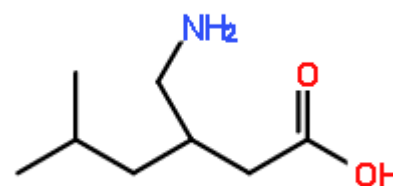

CHROMATOGRAM

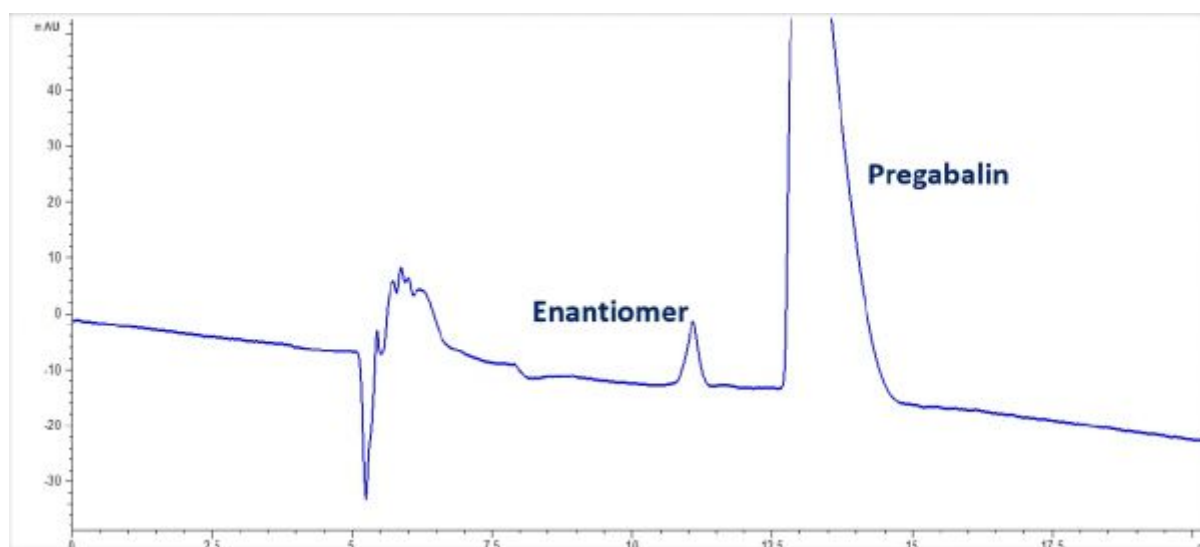

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	4 x 250 mm (3 µm)		
Mobile Phase Composition	water / methanol / acetonitrile / ammonium dihydrogen phosphate (weight) = 100 / 870 / 30 / 0.575		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	10
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 212 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.1		1 - 2	2.37
2	13.0			