


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Rivaroxaban

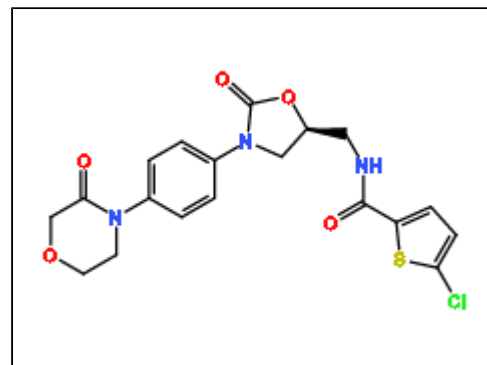

CHROMATOGRAM

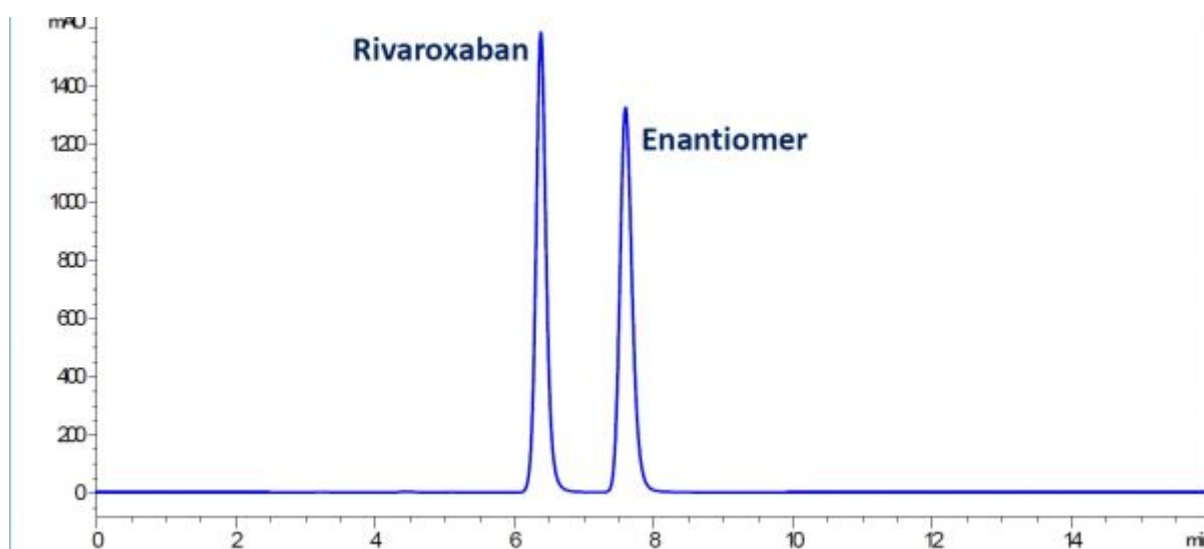

ANALYTICAL REPORT

| | | | |
|-----------------------------|---|----------------------|----|
| Column Trade Name | CHIRALPAK® IH | | |
| Column Size (Particle Size) | 4.6 x 250 mm (5 µm) | | |
| Mobile Phase Composition | ethanol / acetonitrile / diethylamine = 90 / 10 / 0.1 | | |
| Chromatographic Mode | Polar Phase | | |
| Flow rate (mL/min) | 0.7 | Temperature (°C) | 40 |
| Injection Volume (µL) | 0.01 (mg) | Sample Concentration | |
| Detection | UV-VIS 250 nm | | |
| Elution Order | | | |

| No. | Rt | K' | Alpha Value | Resolution |
|-----|-----|----|-------------|------------|
| 1 | 6.4 | | 1 - 2 | 4.20 |
| 2 | 7.6 | | | |