


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ruxolitinib

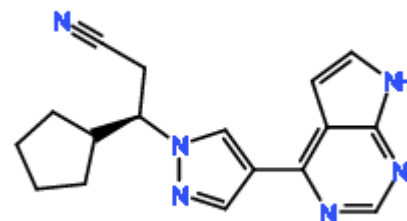

## CHROMATOGRAM

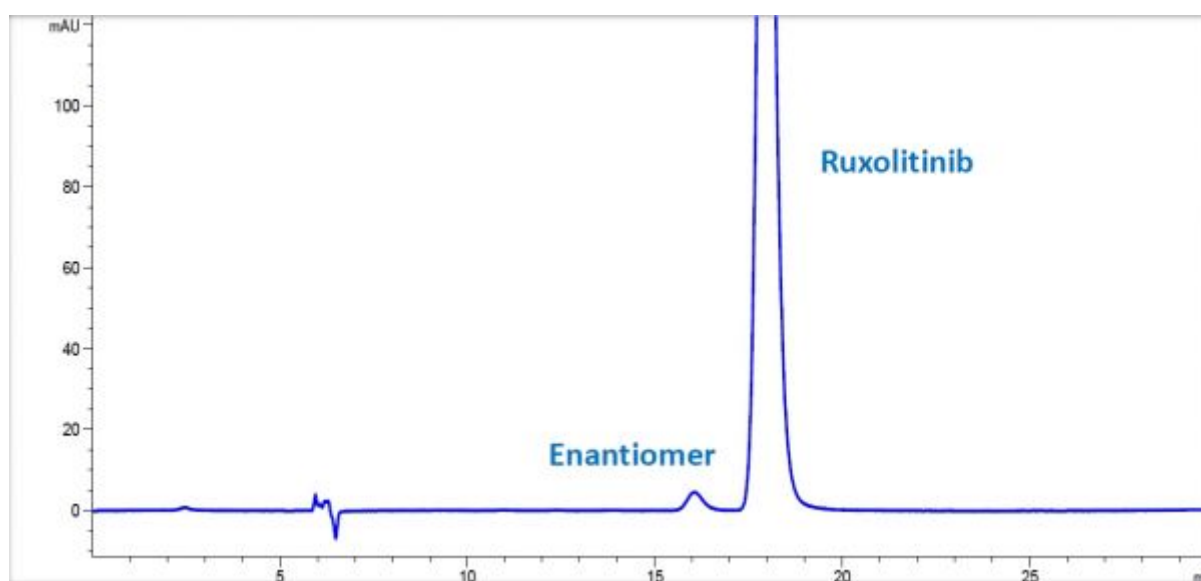

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	ethanol = 100		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	15
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 310 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	16.0		1 - 2	2.37
2	17.9			