


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Salbutamol

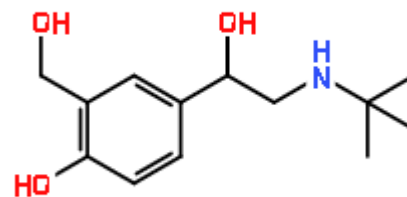

CHROMATOGRAM

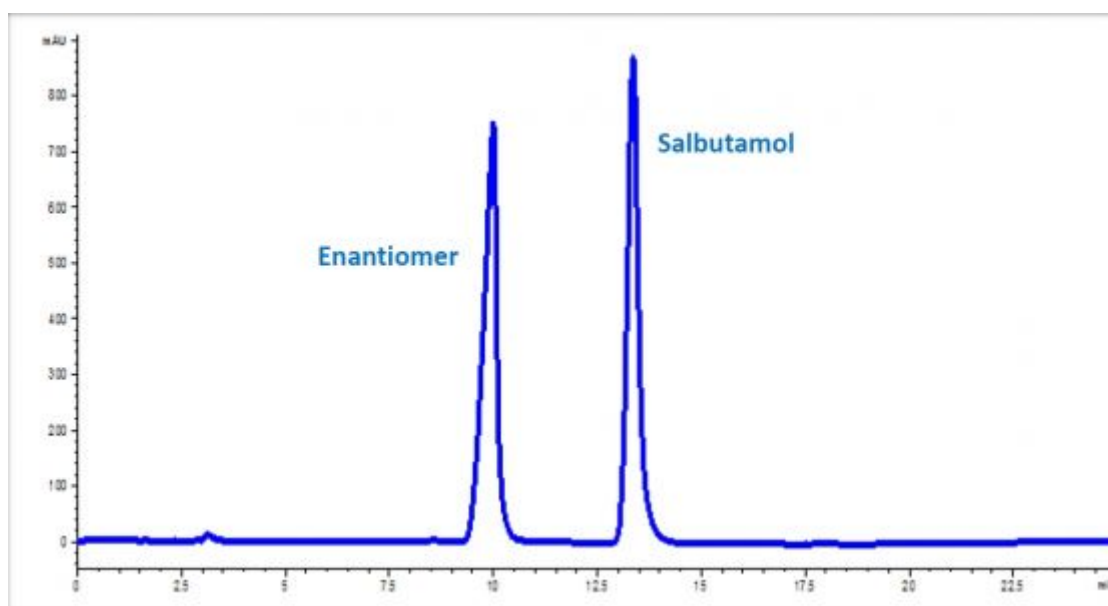

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(-)		
Column Size (Particle Size)	4 x 150 mm (3 µm)		
Mobile Phase Composition	methanol / tetrahydrofuran / water / 2-amino-1-butanol = 490 / 490 / 20 / 1.2		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.0		1 - 2	6.41
2	13.4			