


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Salmeterol

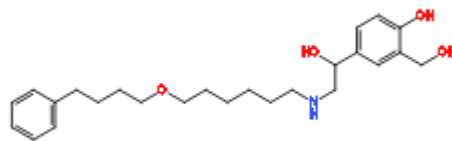

CHROMATOGRAM

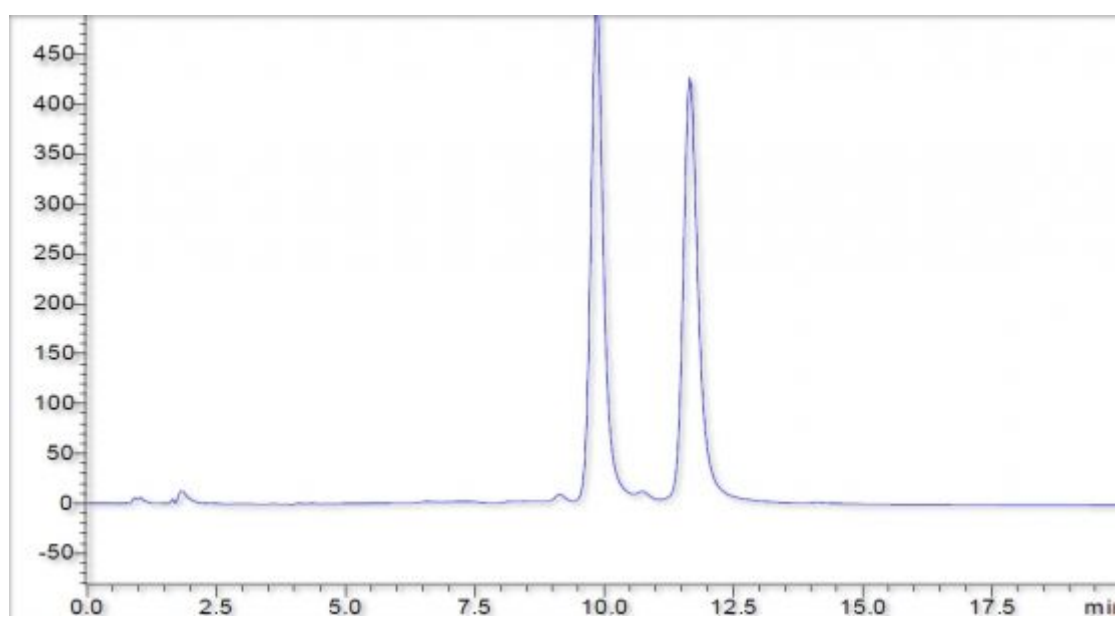

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 3 x 150 mm (3 µm)		
Mobile Phase Composition	methanol / acetonitrile / water / 2-amino-1-butanol = 490 / 490 / 20 / 1.2		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.8		1 - 2	3.43
2	11.6			