


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Stiripentol

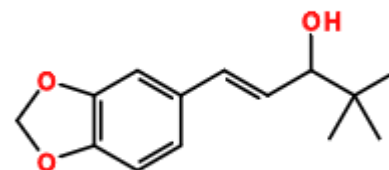

CHROMATOGRAM

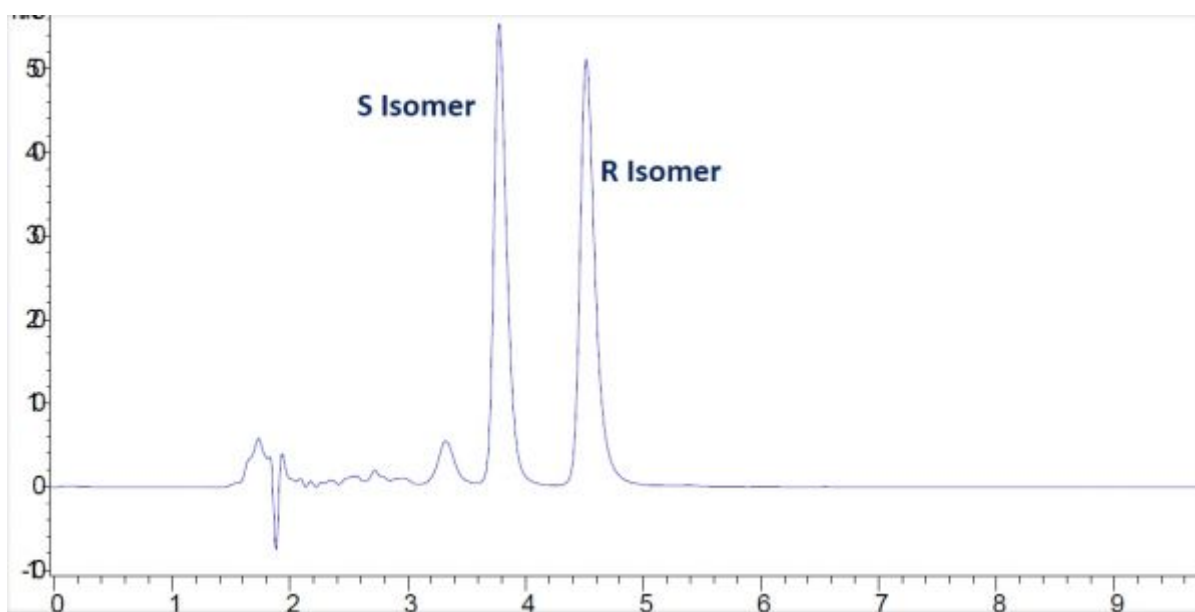

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	10mM NH ₄ OAc aq. / acetonitrile = 30 / 70		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	1 mg/mL
Detection	UV-VIS 268 nm		
Elution Order	1st peak : S, 2nd peak : R		

No.	Rt	K'	Alpha Value	Resolution
1	3.8	1.10	1 - 2	3.57
2	4.5	1.51		