


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Tyrosine

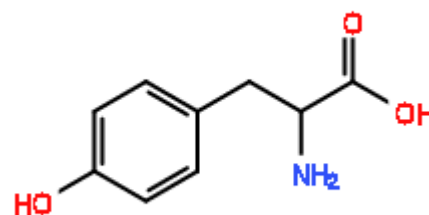

CHROMATOGRAM

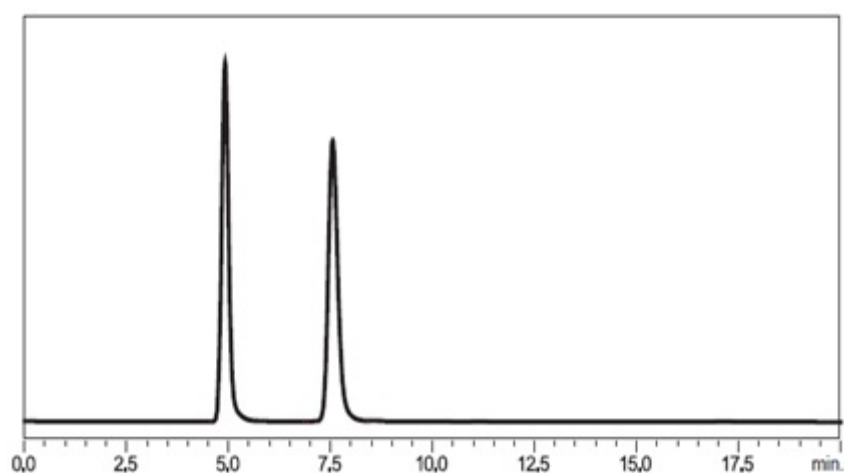

ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.5 / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.2	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 200 nm		
Elution Order			

No.

Rt

K'

Alpha Value

Resolution

1

1 - 2