


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

DL-Asparagine

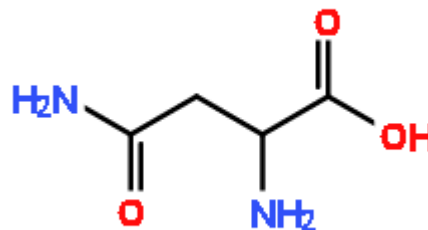

## CHROMATOGRAM

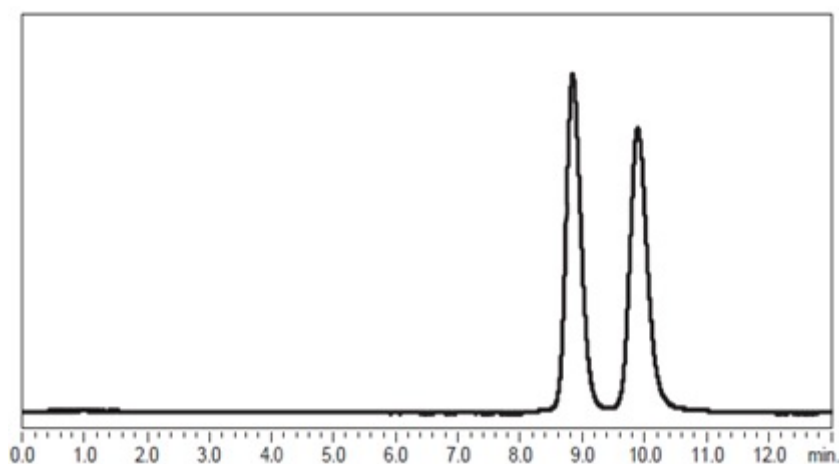

## ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(+)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / methanol = 70 / 30		
Chromatographic Mode	Reverse Phase		

Flow rate (mL/min)	0.1	Temperature (°C)	25
Injection Volume (µL)	0.00025 (mg)	Sample Concentration	
Detection	UV-VIS 200 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1			1 - 2	