


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

1-(1-Naphthyl)-ethanol

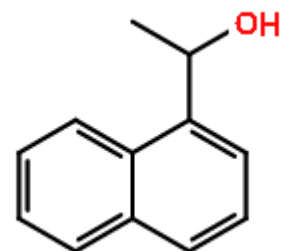

## CHROMATOGRAM

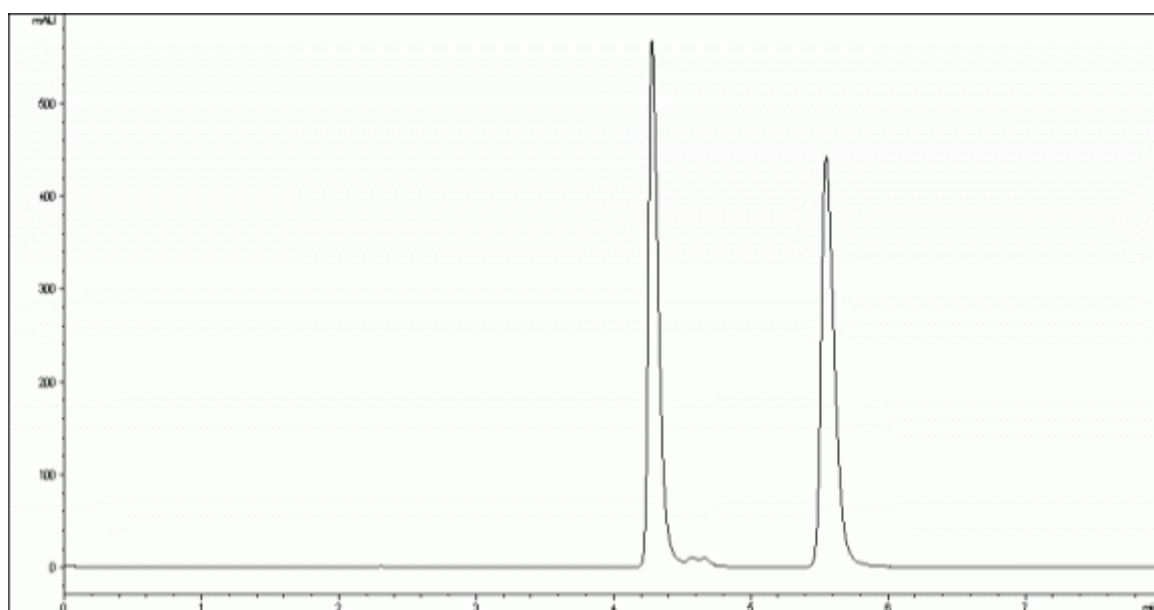

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.3	1.36	1 - 2	8.51
2	5.6	2.06		