


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

2-Aminoadipic acid

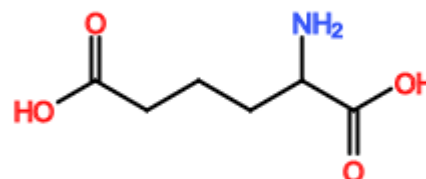

## CHROMATOGRAM

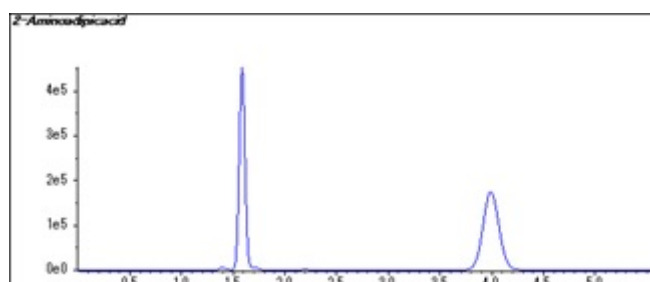

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CROWNPAK® CR-I(+) 3 x 150 mm (5 µm)		
Mobile Phase Composition	acetonitrile / ethanol / water / trifluoroacetic acid = 80 / 15 / 5 / 0.5		
Chromatographic Mode	HILIC		
Flow rate (mL/min)	0.4	Temperature (°C)	30
Injection Volume (µL)		Sample Concentration	
Detection	LC-MS		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.6		1 - 2	13.03
2	4.0			

## Comment

This data was provided by Professor Fukusaki of Osaka University.

References ;

Novel high-throughput and widely-targeted liquid chromatography–time of flight mass spectrometry method for d-amino acids in foods

Yutaka Konya, Moyu Taniguchi, Eiichiro Fukusaki

Volume 123, Issue 1, January 2017, Pages 126-133