


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Guaifenesin

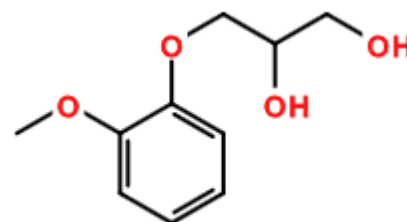

CHROMATOGRAM

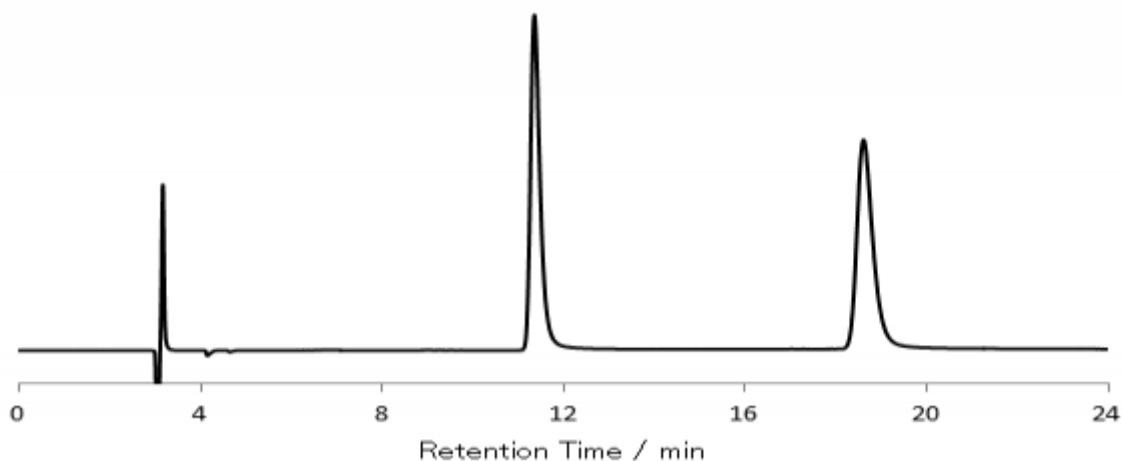

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB N-5		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / tetrahydrofuran = 70 / 30		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.4	2.90	1 - 2	14.56
2	18.6	5.30		