


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Pindolol

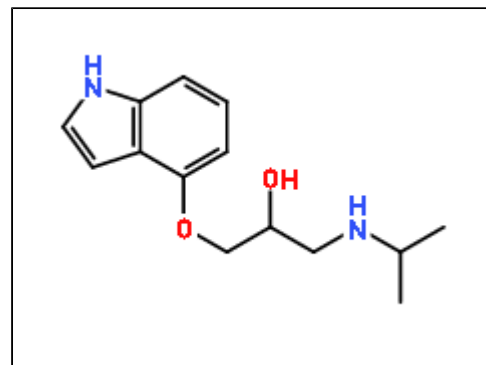

CHROMATOGRAM

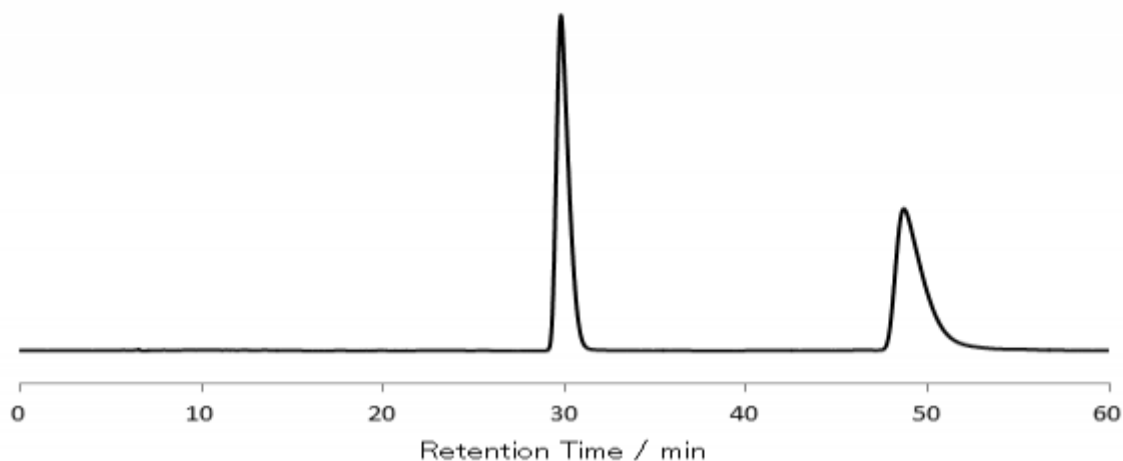

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB N-5		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	0.1M Potassium hexafluorophosphate aq. / acetonitrile = 70 / 30		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	29.8	4.00	1 - 2	9.64
2	48.7	7.20		