

APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Dimetotiazine



CHROMATOGRAM

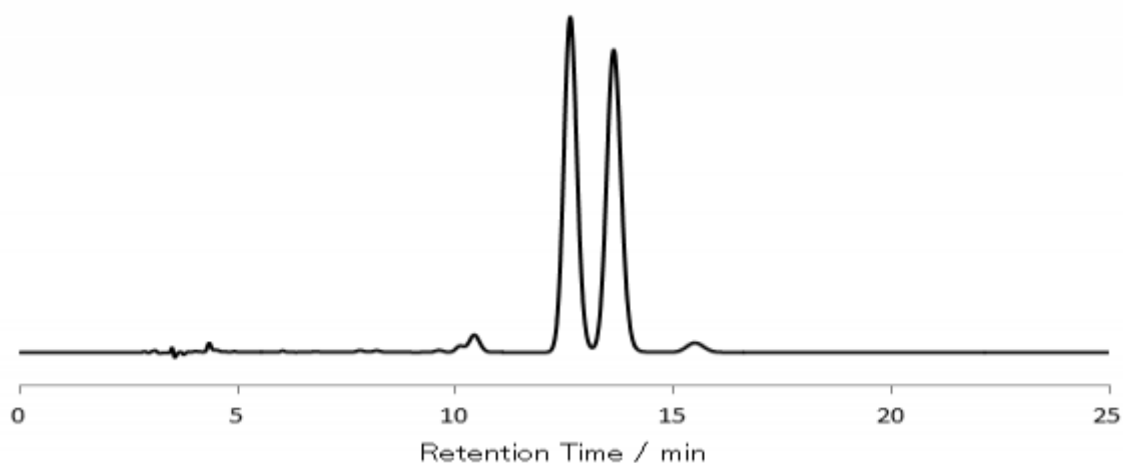

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (μL)		Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.6	3.20	1 - 2	1.62
2	13.6	3.50		