


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

5-Methyl-5-phenylhydantoin

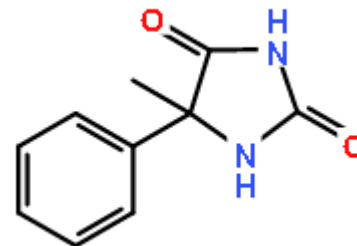

CHROMATOGRAM

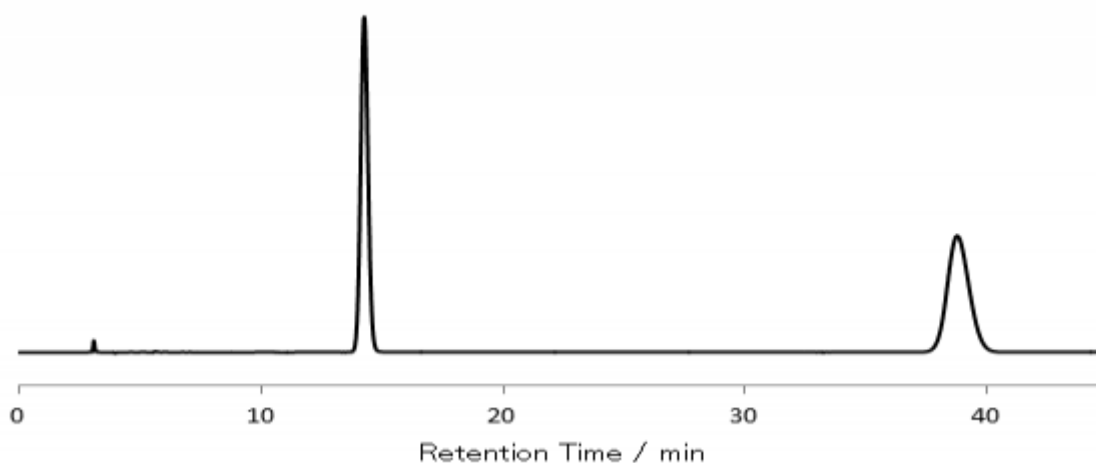

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μ m)		
Mobile Phase Composition	n-hexane / ethanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature ($^{\circ}$ C)	25
Injection Volume (μ L)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	14.3	3.80	1 - 2	21.62
2	38.8	11.90		