


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ofloxacin

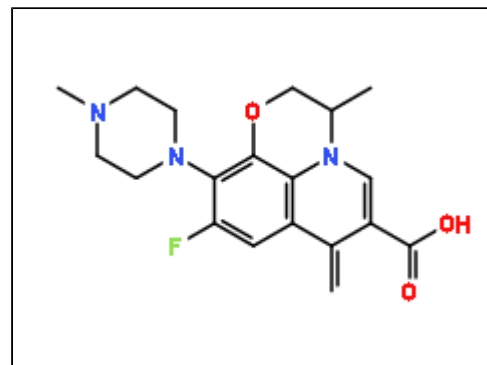

CHROMATOGRAM

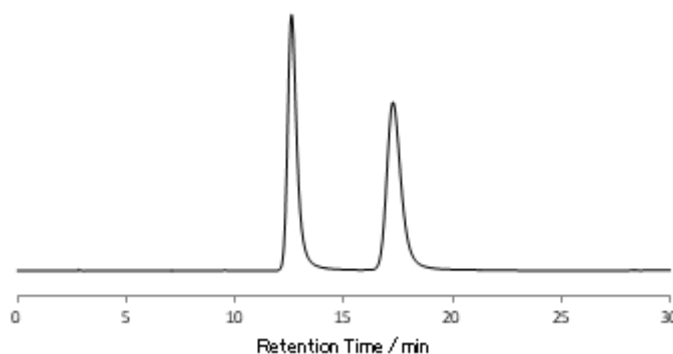

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IH 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / acetic acid / ethylenediamine = 70 / 30 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 290 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.6	3.20	1 - 2	4.83
2	17.2	4.70		