


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Valine

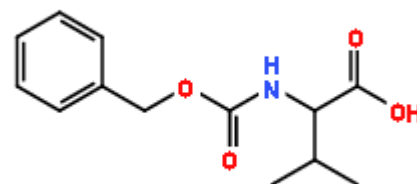

CHROMATOGRAM

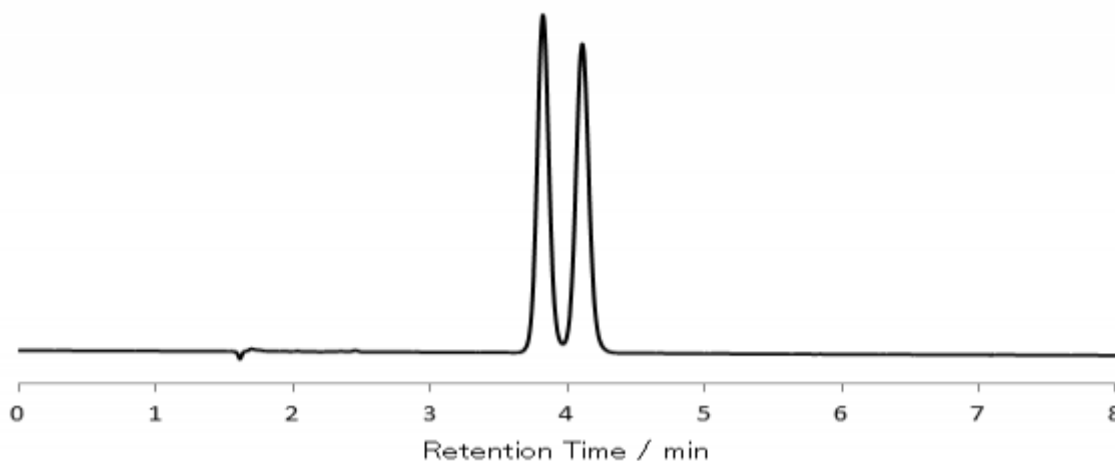

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	formic acid aq. pH2.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.003 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.8	1.27	1 - 2	1.70
2	4.1	1.44		