


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

5-Methyl-5-phenylhydantoin

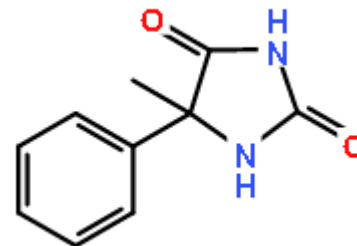

## CHROMATOGRAM

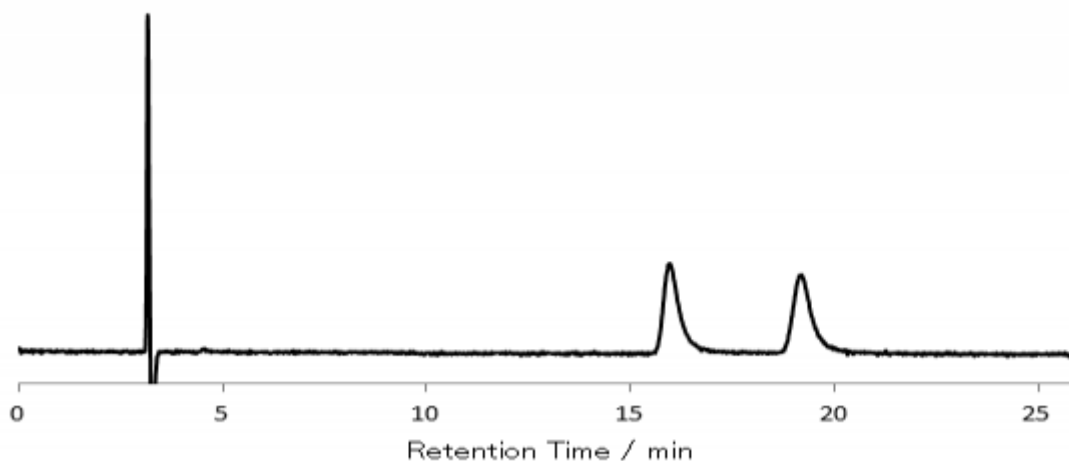

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / tetrahydrofuran = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	16.0	4.32	1 - 2	4.53
2	19.2	5.40		