


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Flurbiprofen

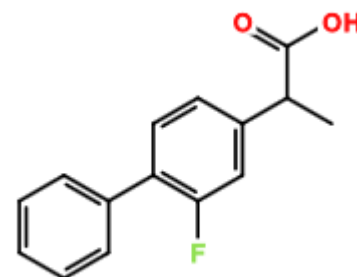

CHROMATOGRAM

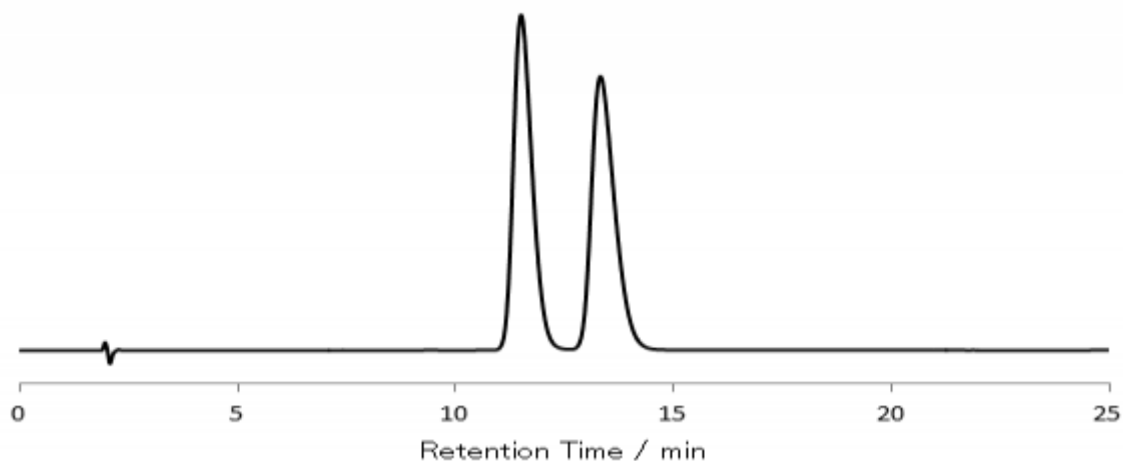

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	formic acid aq. pH2.0 / methanol = 20 / 80		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 245 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	11.5	4.75	1 - 2	2.07
2	13.3	5.66		