


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Hydroxyzine

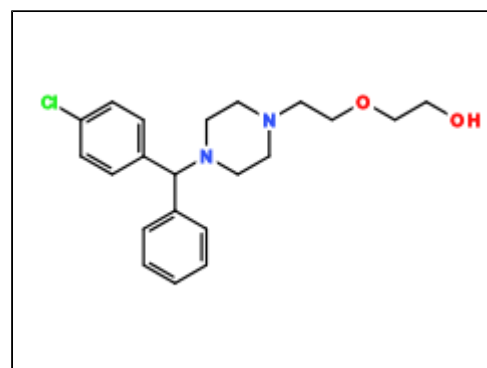

CHROMATOGRAM

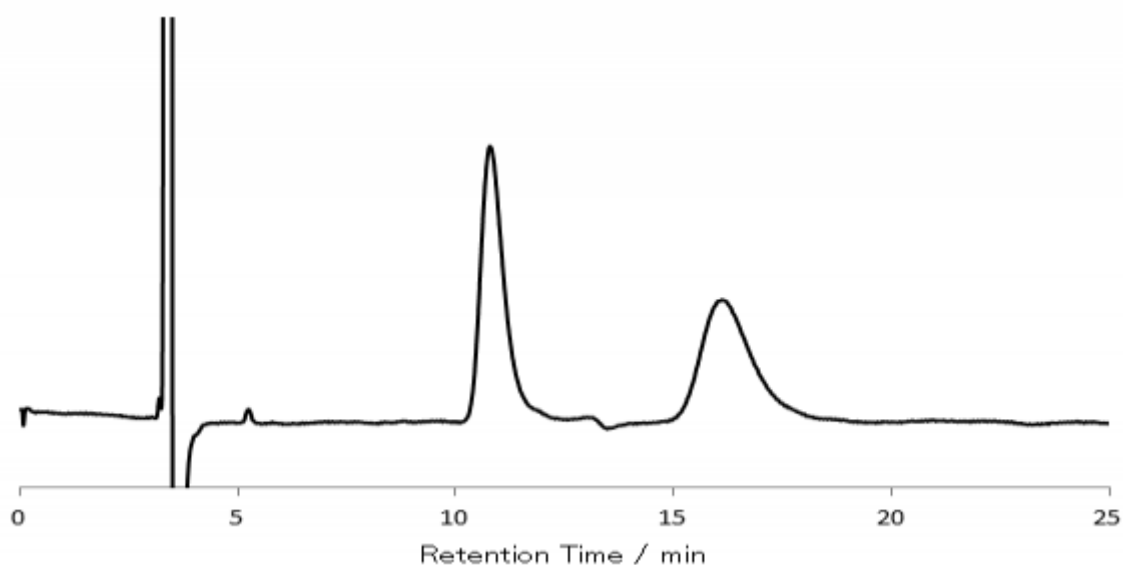

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / tetrahydrofuran / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 236 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.8	2.59	1 - 2	3.48
2	16.1	4.37		