


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

5-Methyl-5-phenylhydantoin

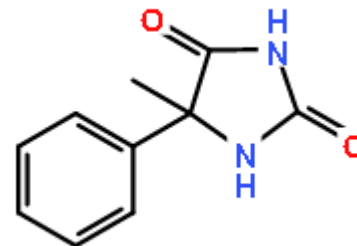

CHROMATOGRAM

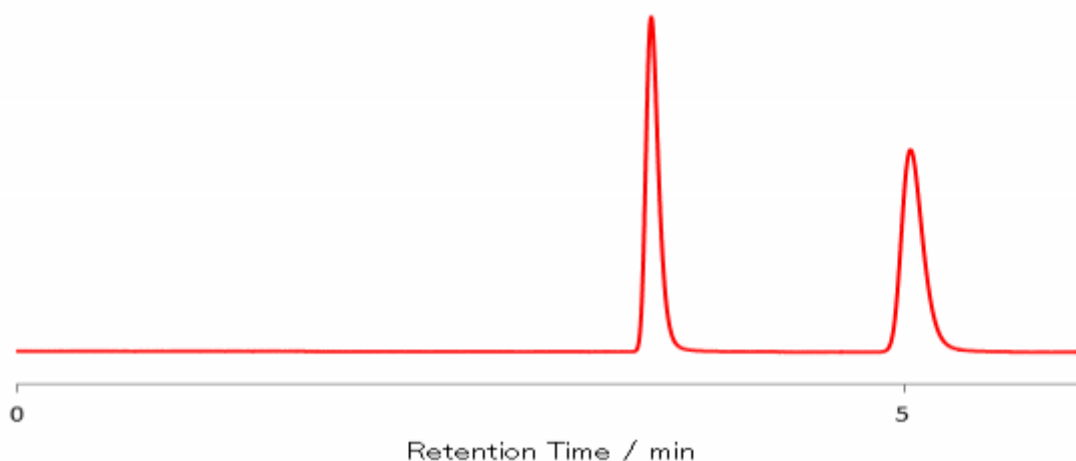

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	water / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.21	Temperature (°C)	25
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 205 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.6	0.68	1 - 2	8.17
2	5.0	1.37		