


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Naproxen

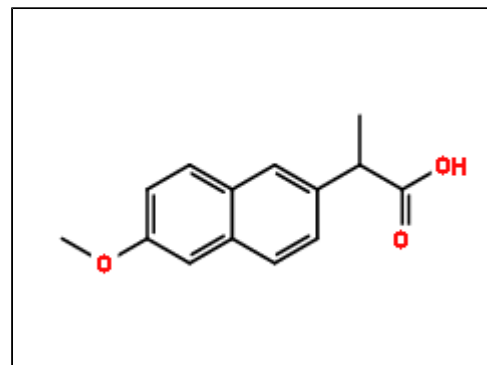

CHROMATOGRAM

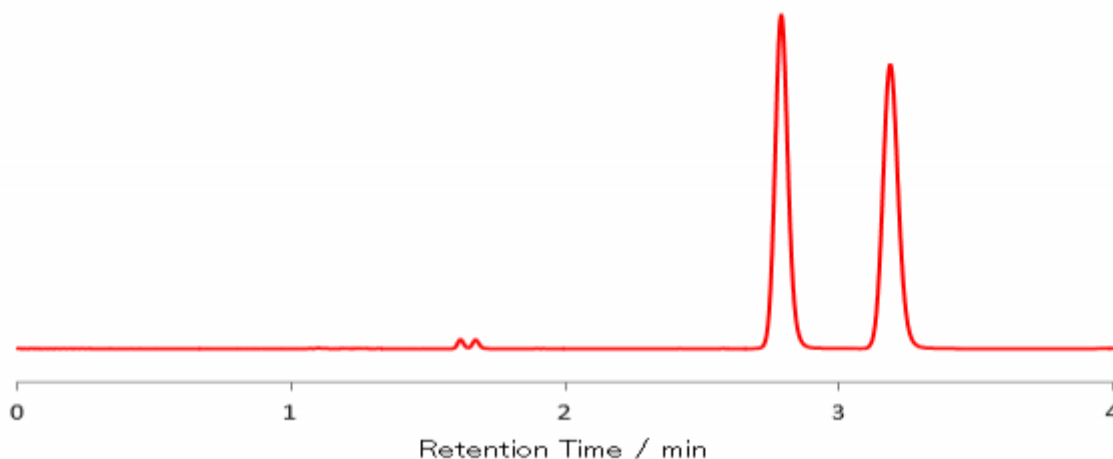

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.425	Temperature (°C)	25
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.8	1.60	1 - 2	4.15
2	3.2	1.97		