


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1,1'-Bi-2-naphthol

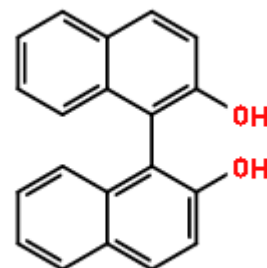

CHROMATOGRAM

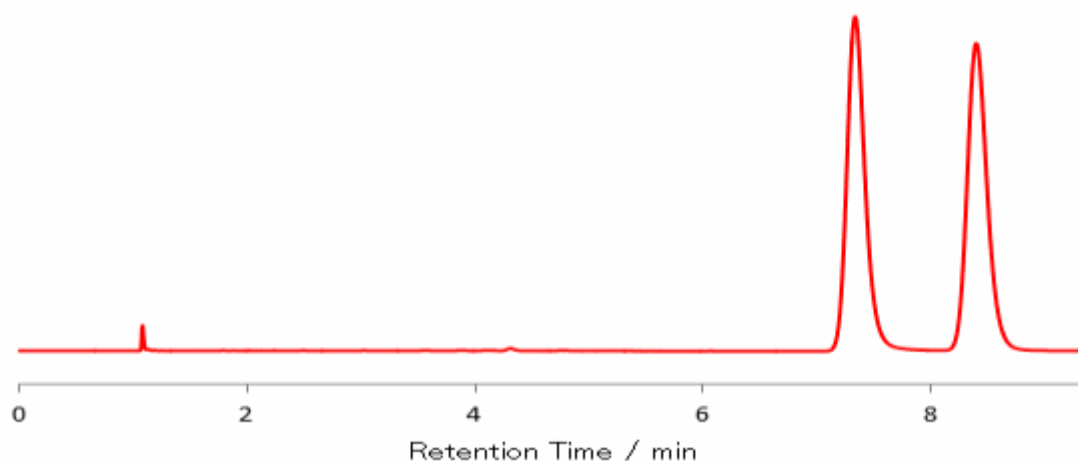

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	n-hexane / ethanol = 95 / 5		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.4	Temperature (°C)	25
Injection Volume (µL)	0.0001 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.3	5.79	1 - 2 1.17	3.55
2	8.4	6.77		