


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Hexobarbital

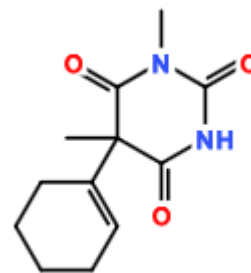

CHROMATOGRAM

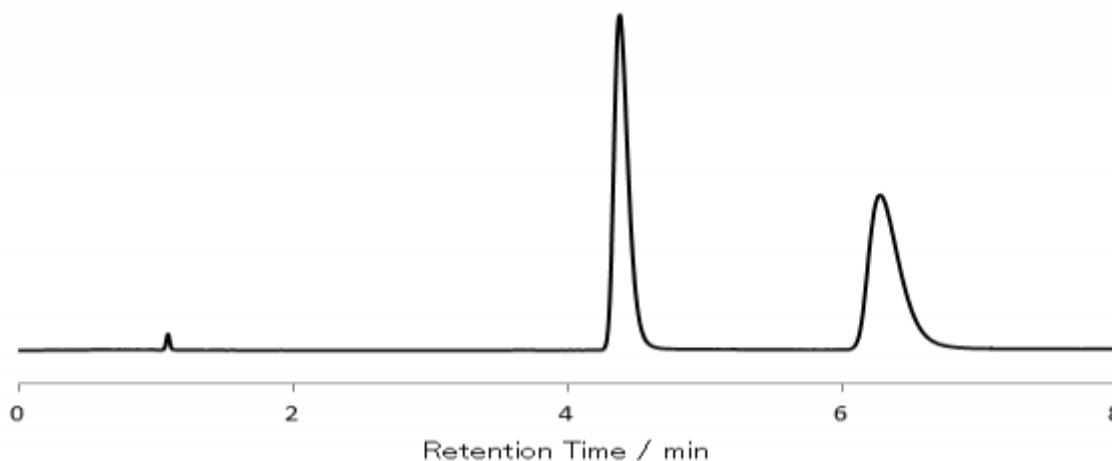

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.425	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 205 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.4	3.24	1 - 2	6.23
2	6.3	5.08		