


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cobalt(III) tris(acetylacetonate)

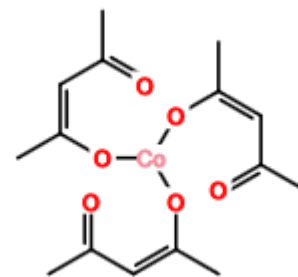

## CHROMATOGRAM

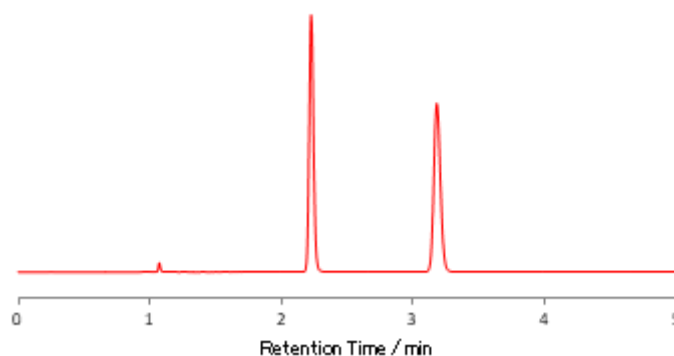

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK <sup>®</sup> IC-U		
Column Size (Particle Size)	3 x 100 mm (1.6 μm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	0.425	Temperature (°C)	25
Injection Volume (μL)	0.00004 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.2	1.08	1 - 2	12.53
2	3.2	1.97		