


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Prilocaine

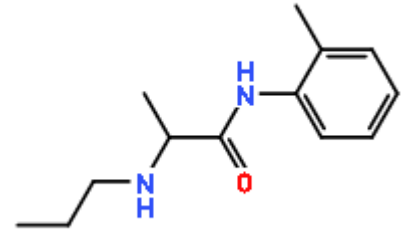

CHROMATOGRAM

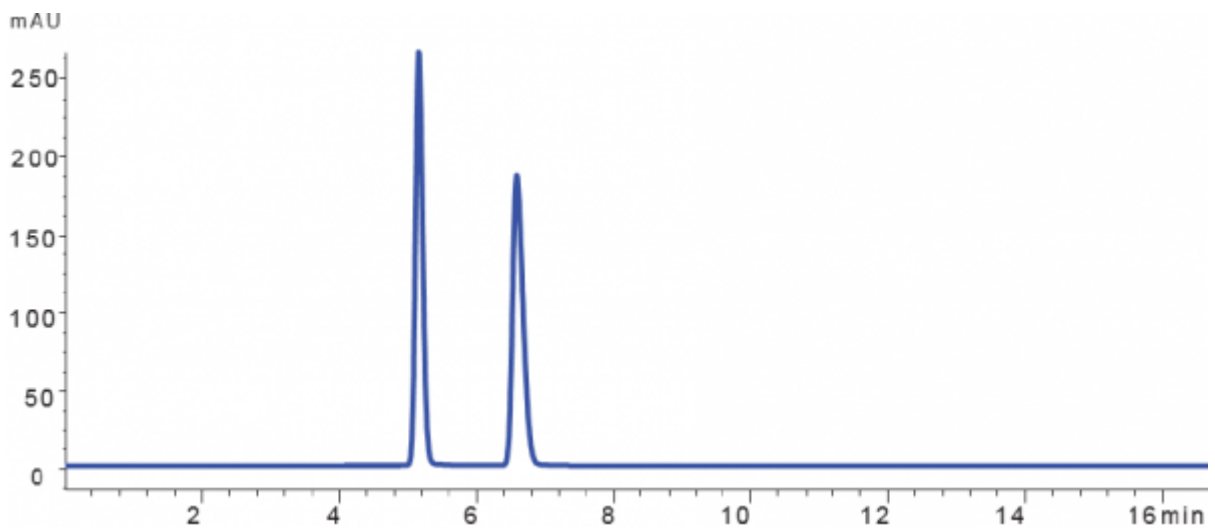

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IG-3 4.6 x 150 mm (3 µm)		
Mobile Phase Composition	100mM ammonium bicarbonate aq. / acetonitrile = 5 / 95		
Chromatographic Mode	HILIC		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.3		1 - 2	6.10
2	5.5			