


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Elagolix

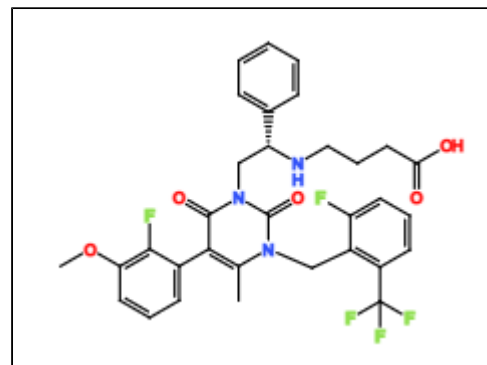

CHROMATOGRAM

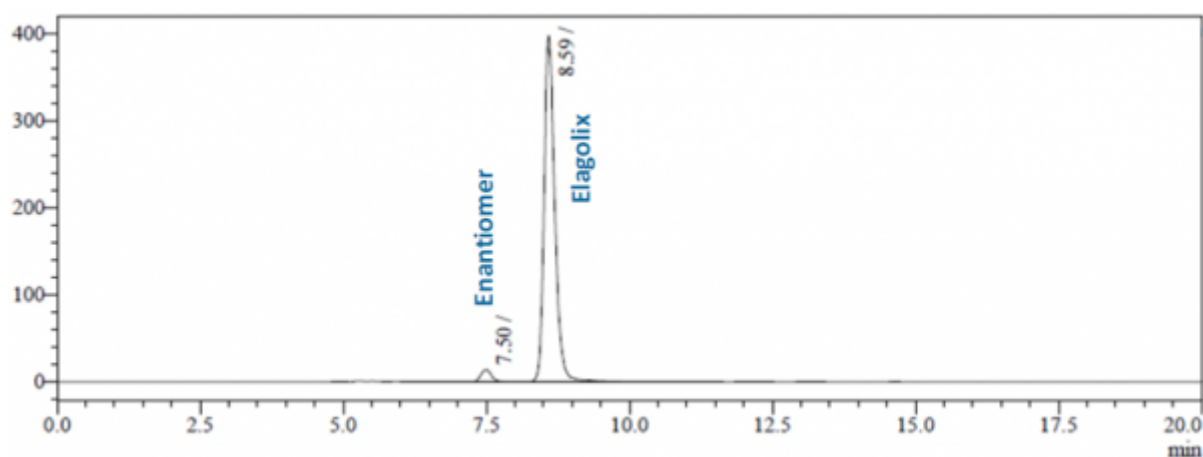

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 4 x 250 mm (3 µm)		
Mobile Phase Composition	methanol / water / acetic acid / diethylamine = 980 / 20 / 2.9 / 2.6		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	40
Injection Volume (µL)	10	Sample Concentration	0.5 mg/mL in MeOH
Detection	UV-VIS 275 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.5		1 - 2	3.20
2	8.6			