


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Molnupiravir

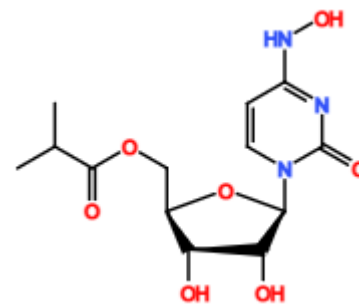

CHROMATOGRAM

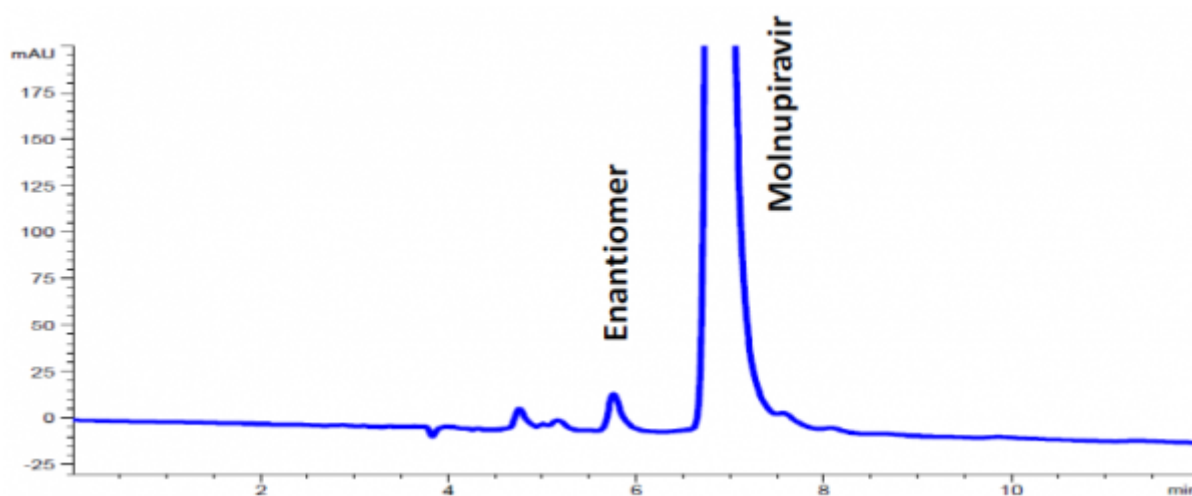

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IC-3 4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / ethanolamine = 50 / 25 / 25 / 0.3		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	0.5 mg/mL in Mobile phase
Detection			
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.8		1 - 2	4.10
2	6.8			