


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Troxiptide

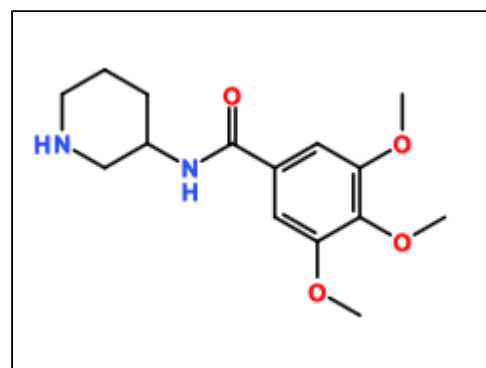

## CHROMATOGRAM

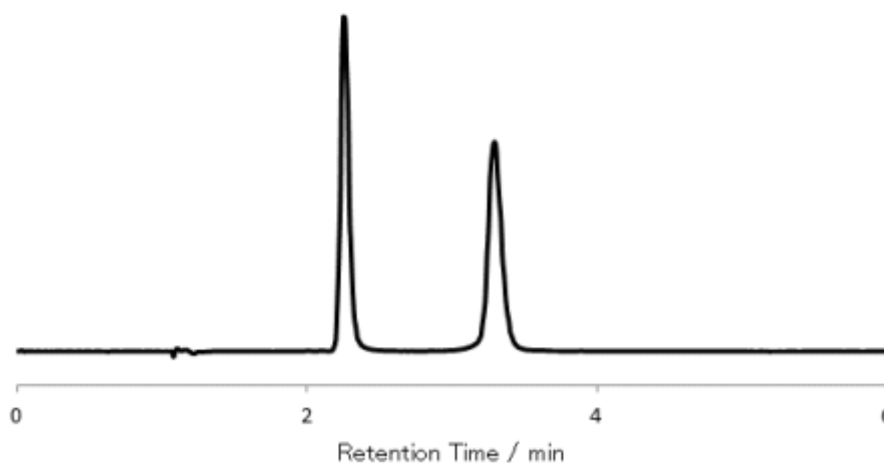

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO <sub>2</sub> / methanol / water / diethylamine = 60 / 40 / 1 / 0.1		
Chromatographic Mode	SFC		
Flow rate (mL/min)	2.0	Temperature (°C)	40
Injection Volume (µL)	5.0	Sample Concentration	1.0 mg/mL in MeOH
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.3	0.50	1 - 2	2.37
2	3.3	1.19		7.58