


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ethotoin

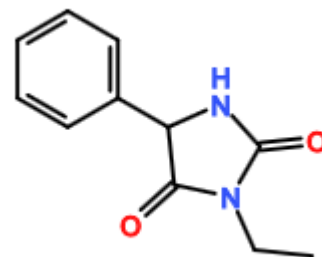

CHROMATOGRAM

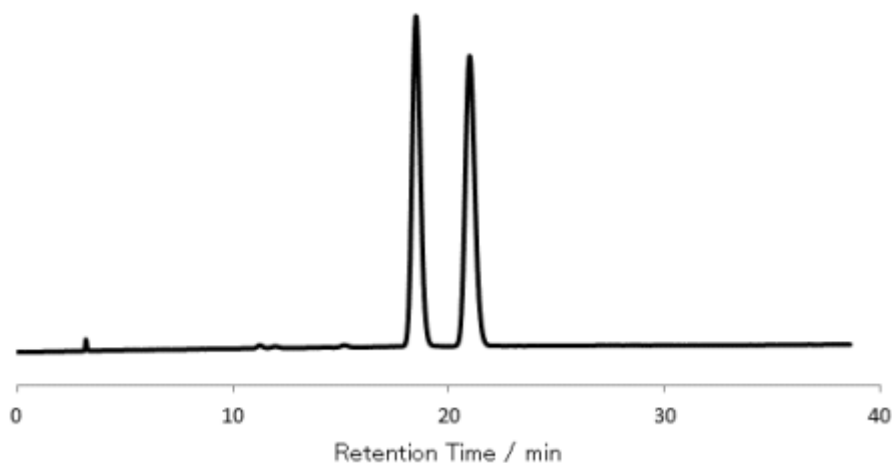

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5.0	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 200 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	18.5	5.15	1 - 2	3.20
2	20.9	5.98		